


**Organochlorine Pesticides (EPA 8081) Mix # 1 on ZB-35**

**Column:** Zebtron™ ZB-35, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** 35% Phenyl 65% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G003-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min  
**Carrier Gas:** Constant Flow Helium, 0.8 mL/min  
**Injection:** On-Column :1 0.5 µL @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes 10 ppm in Acetone

**Products used in this application:**

18394

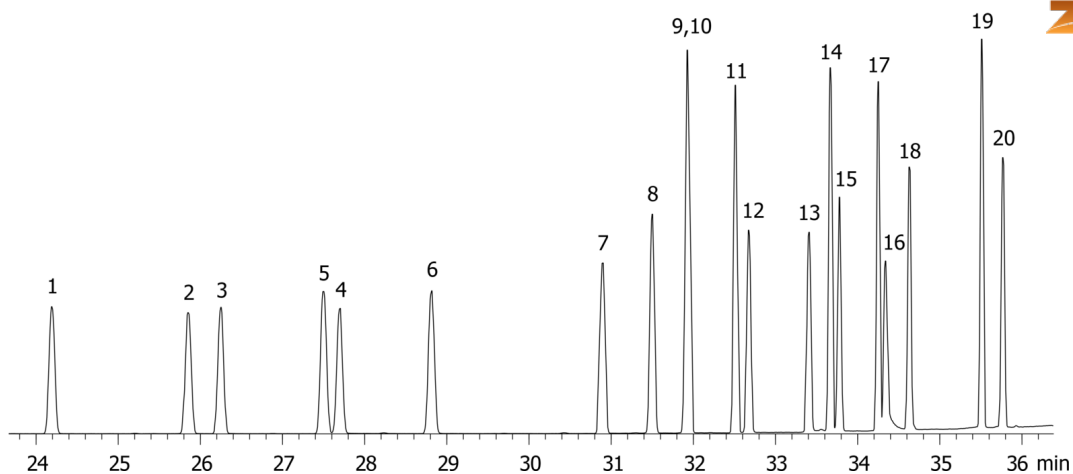

## Organochlorine Pesticides (EPA 8081) Mix # 1 on ZB-35

### ANALYTES:

- 1 alpha-BHC
- 2 gama-BHC
- 3 b-BHC
- 4 delta-BHC
- 5 Heptachlor
- 6 Aldrin
- 7 Heptachlor epoxide
- 8 gamma-Chlordane (trans)
- 9 alpha-Chlordane (cis)
- 10 Endosulfan I
- 11 DDE
- 12 Dieldrin
- 13 Endrin
- 14 DDD
- 15 Endosulfan II
- 16 Endrin aldehyde
- 17 DDT
- 18 Endosulfan sulfate
- 19 Methoxychlor
- 20 Endrin ketone

