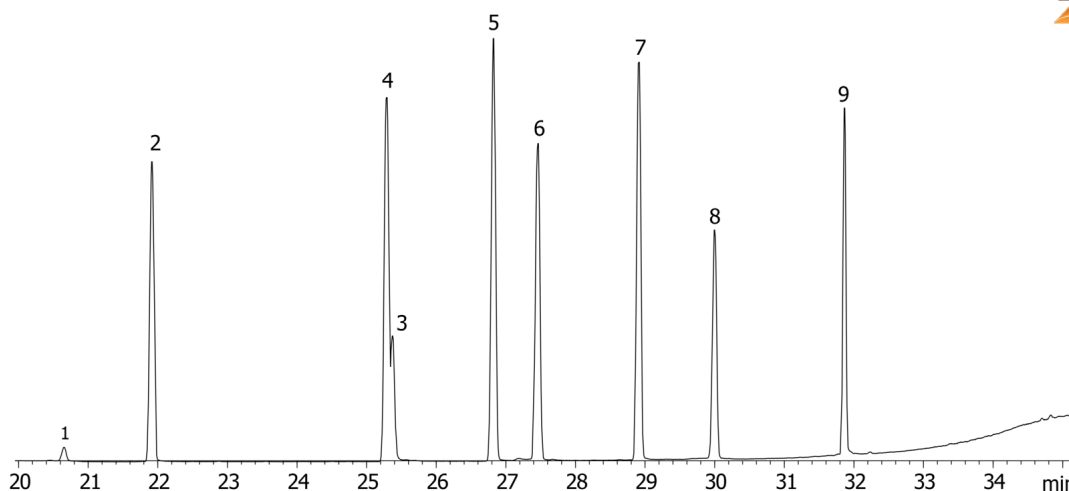


Organophosphorous Pesticides (EPA 8141) Mix # 4 on ZB-1701P

Column: Zebron[™] ZB-1701, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 14% Cyanopropylphenyl 86% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G006-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 290 °C @ 15 °C/min for 25 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: On-Column :1 0.5 µL @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes 10 ppm in Acetone

**Products used in this application:**

18435

**ANALYTES:**

- 1 Trichlorfon
- 2 Thionazin
- 3 Dicrotophos
- 4 Fonofos
- 5 Dichlofenthion
- 6 Chlorpyrifos-methyl
- 7 Aspon
- 8 Fenitrothion
- 9 Crotoxyphos

