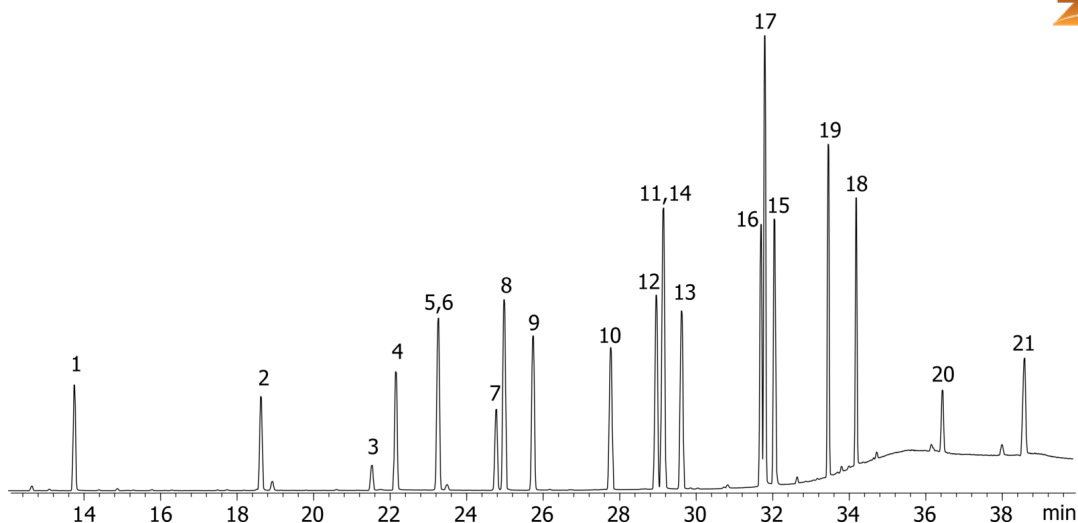


Organophosphorous Pesticides (EPA 8141) Mix # 1 on ZB-1701P

Column: Zebtron™ ZB-1701, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 14% Cyanopropylphenyl 86% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G006-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 290 °C @ 15 °C/min for 25 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: On-Column :1 0.5 µL @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes 10 ppm in Acetone

**Products used in this application:**

Organophosphorous Pesticides (EPA 8141) Mix # 1 on ZB-1701P

ANALYTES:

- 1 Dichlorvos
- 2 Mevinphos
- 3 Demeton
- 4 Ethoprop
- 5 Naled
- 6 Phorate
- 7 Demeton
- 8 Diazinon
- 9 Disulfoton
- 10 Fenchlorphos
- 11 Methyl parathion
- 12 Chlorpyrifos
- 13 Fenthion
- 14 Trichloronate
- 15 Tetrachlorvinphos
- 16 Prothiophos
- 17 Merphos
- 18 Fensulfothion
- 19 Sulprofos
- 20 Azinphos-methyl
- 21 Coumaphos
- 22 Tetrachlorvinphos (stirofos)
- 23 Tokuthion (prothiophos)
- 24 Merphos
- 25 Fensulfothion
- 26 Bolstar (sulprofos)
- 27 Azinphos-methyl (guthion)
- 28 Coumaphos

