


Zebron ZB-Drug-1 Test Mix (AG0-8431)**Column:** ZebronTM ZB-Drug-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G023-11**Oven Profile:** 140 °C Isothermal**Carrier Gas:** Constant Flow Hydrogen, 1.12 mL/min**Injection:** Split 89:1 1 µL @ 250°C**Detection:** Flame Ionization (FID) (325°C)**Analyst Note:** Analytes are 250ug/mL in Acetone (PN: AG0-8431)**Products used in this application:**

18461

**ANALYTES:**

- 1 Dodecane
- 2 Tridecane
- 3 4-Chlorophenol
- 4 Tetradecane
- 5 1-Undecanol
- 6 1-Methylnaphthalene
- 7 Dicyclohexylamine
- 8 Hexadecane (C16)

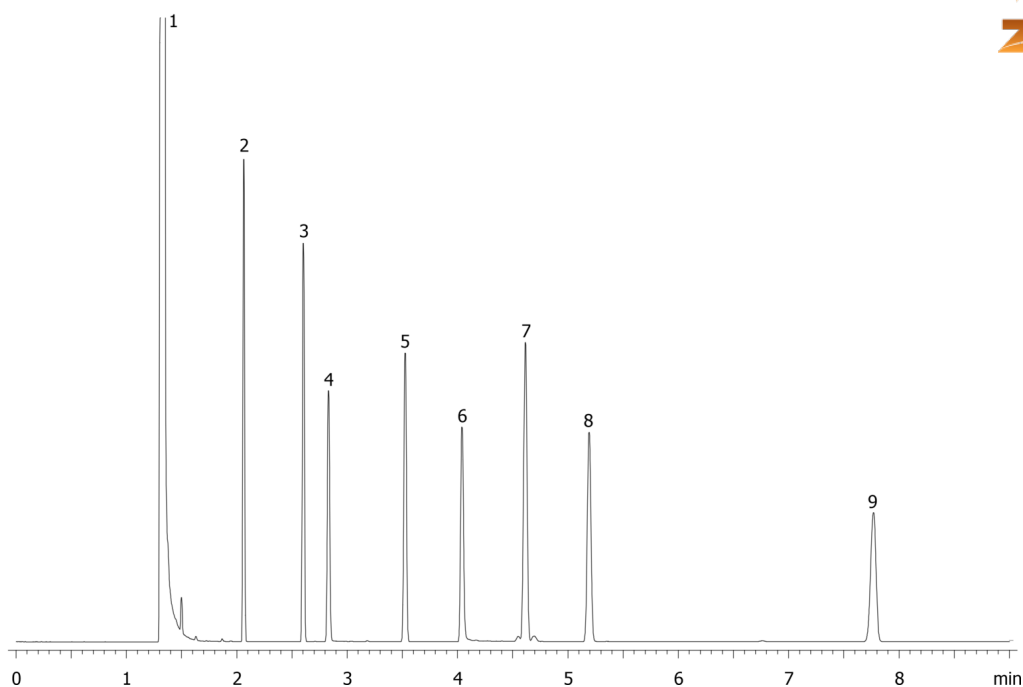