

Malachite Green and Related Dyes on Luna HST 2.5um C18(2) 50x2.0 Column

Column: Luna® 2.5 µm C18(2) 100 Å, LC Column 50 x 2 mm, Ea

Dimensions: 50 x 2 mm ID

Order No: 00B-4446-B0

Elution Type: Isocratic

Eluent A: 0.1% Formic acid in water

Eluent B: 0.1% Formic acid in 75:25 acetonitrile:methanol


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	5	55	45	

Flow Rate: 0.4 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

18464

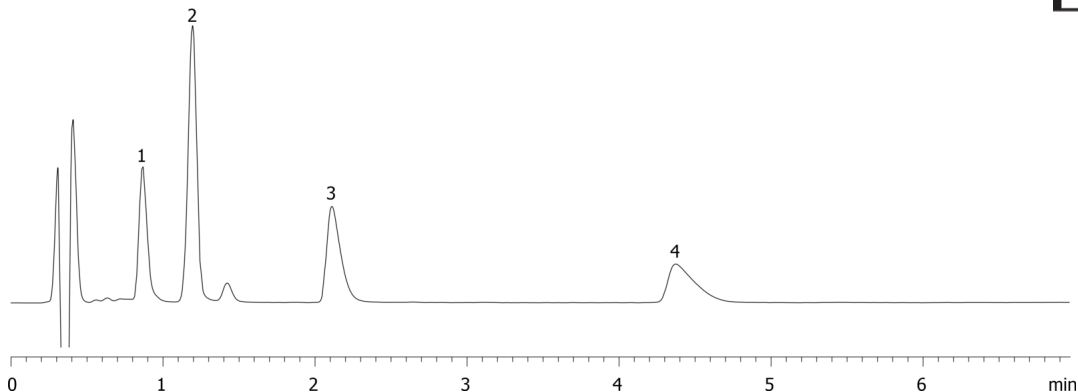

Products used in this application:



ANALYTES:

- 1** Malachite Green
Retention Time: 0.866 min
- 2** Leucomalachite Green
Retention Time: 1.194 min
- 3** Crystal violet
Retention Time: 2.11 min
- 4** Brilliant Green
Retention Time: 4.37 min

