


BioMethanol Analysis - PME on ZB-FFAP**Column:** Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G009-11**Oven Profile:** 50 °C for 8 min**Carrier Gas:** Constant Flow Hydrogen, 1 mL/min**Injection:** Split 50:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Courtesy of ASG (www.asg-analytik.de), a specialized Biofuel Laboratory!**Products used in this application:****ANALYTES:**

- 1 Hexaninoel
- 2 Methanol-FAME

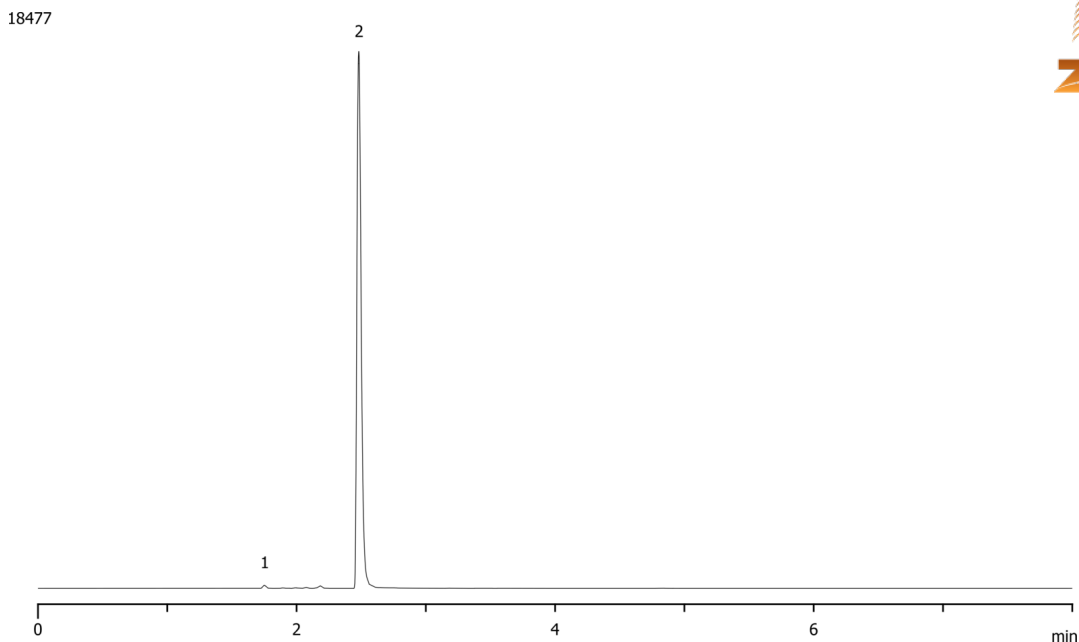