

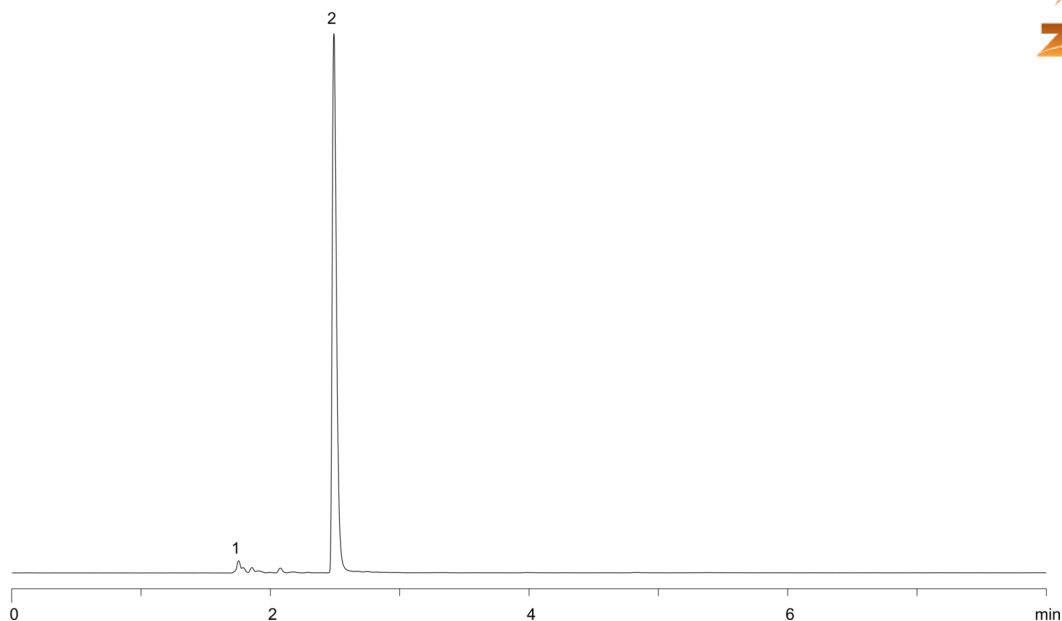
BioMethanol Analysis - SME on ZB-FFAP

Column: Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: Nitroterephthalic Acid Modified Polyethylene Glycol
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G009-11
Oven Profile: 50 °C for 8 min

Carrier Gas: Constant Flow Hydrogen, 1 mL/min
Injection: Split 50:1 1 μ L @ 250°C
Detection: Flame Ionization (FID) (300°C)
Analyst Note: Courtesy of ASG (www.asg-analytik.de), a specialized Biofuel Laboratory!

**Products used in this application:**

18479

**ANALYTES:**

- 1 Hexaninoel
- 2 Methanol-FAME

