

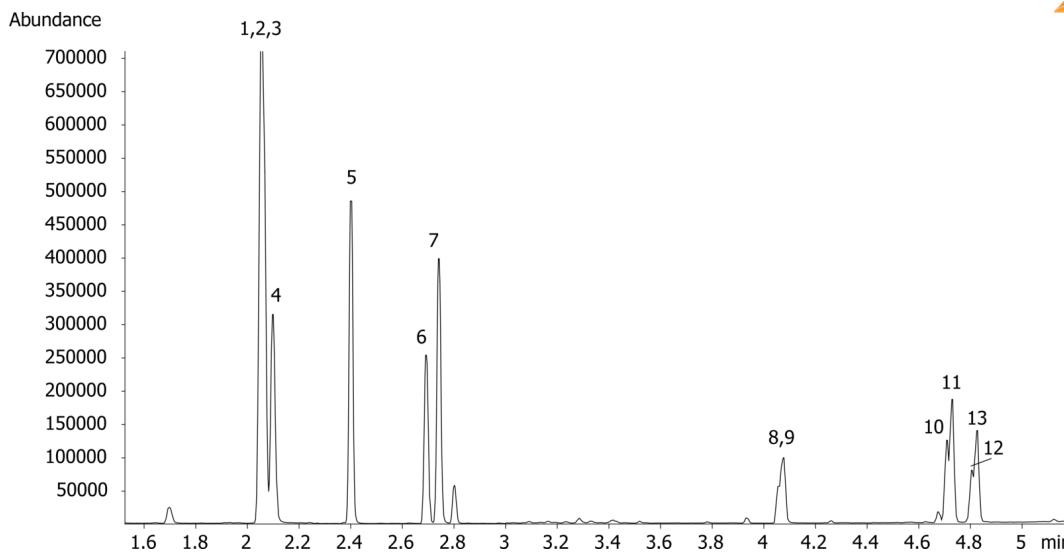
Alternative amphetamines by MSD on a ZB-Drug1**Column:** Zebron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:** 90 °C to 320 °C @ 20 °C/min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10:1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 20ppm derivatized with HFBA

Sim ions:

(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):

Methamphetamine (*254,210,118): D5-MDA (*167,380,244): MDA (*162,375,240): D5-MDMA (*258,213,164):

MDMA (*254,210,389): D5-MDEA (*273,241,408): MDEA (*268,240,403): Phentermine (*254,214,132):

**Products used in this application:**

Alternative amphetamines by MSD on a ZB-Drug1

ANALYTES:

- 1 Phenylpropanolamine
- 2 Phentermine
- 3 Amphetamine-d11
- 4 Amphetamine
- 5 Ephedrine
- 6 Methamphetamine-d14
- 7 Methamphetamine
- 8 MDA-d5
- 9 MDA
- 10 MDMA-d5
- 11 MDMA
- 12 MDEA-d5
- 13 MDEA

