

## Steroids on Luna 3µm PFP(2) 50 x 2mm ID (in acetonitrile)

**Column:** Luna® 3 µm PFP(2) 100 Å, LC Column 50 x 2 mm, Ea

**Dimensions:** 50 x 2 mm ID

**Order No:** 00B-4447-B0

**Elution Type:** Gradient

**Eluent A:** Water

**Eluent B:** CH<sub>3</sub>CN

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 80    | 20    |
|                   | 2        | 10         | 50    | 50    |
|                   | 3        | 10.1       | 80    | 20    |
|                   | 4        | 12         | 80    | 20    |

**Flow Rate:** 0.8 mL/min

**Col. Temp.:** 25 °C


**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (25 °C)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, PFP 4 x 2.0mm ID, 10/Pk Part No.: AJ0-8326
- Holder Part No.: KJ0-4282



Products used in this application:

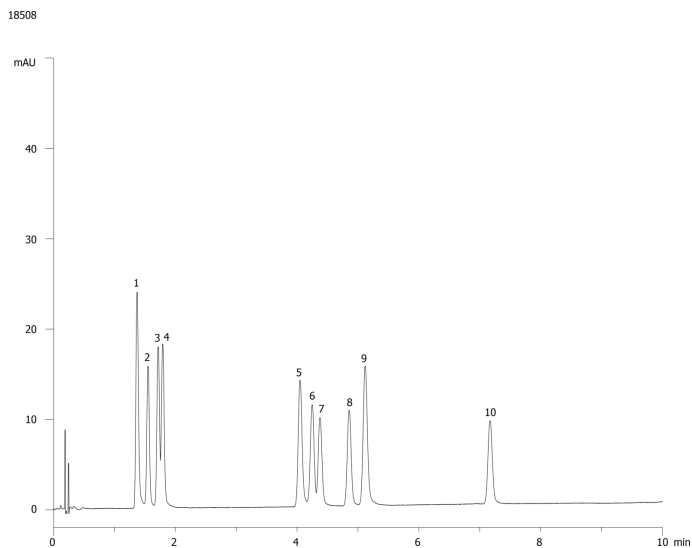

## Steroids on Luna 3 $\mu$ m PFP(2) 50 x 2mm ID (in acetonitrile)

### ANALYTES:

- 1 Estriol
- 2 Hydrocortisone
- 3 Prednisone
- 4 Cortisone
- 5 Estradiol
- 6 Cortisone acetate
- 7 21-Hydroxyprogesterone
- 8 17-Hydroxyprogesterone
- 9 Estrone
- 10 Deoxycorticosterone

