


US PAHs on ZB-35 by MSD

Column: Zebron[™] ZB-35, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: 35% Phenyl 65% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G003-11
Oven Profile: 80 °C for 0.66 to 250 °C @ 20 °C/min to 300 °C @ 8 °C/min to 360 °C @ 20 °C/min for 6 min
Carrier Gas: Constant Flow Helium, 1.2 mL/min
Injection: On-Column :1 μ L @ 83°C
Detection: Mass Selective (MSD) (360°C)
Analyst Note: Analytes are 10 ppm in dichloromethane

**Products used in this application:**

18523

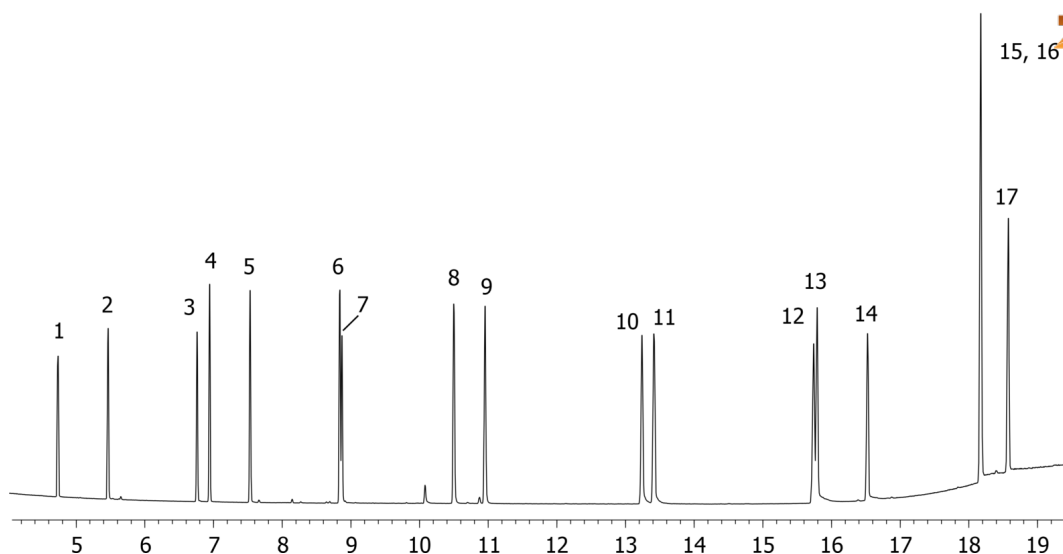

US PAHs on ZB-35 by MSD

ANALYTES:

- 1 Naphthalene
- 2 2-Methylnaphthalene
- 3 Acenaphthylene
- 4 Acenaphthene
- 5 Fluorene
- 6 Phenanthrene
- 7 Anthracene
- 8 Fluoranthene
- 9 Pyrene
- 10 Benz[a]anthracene
- 11 Chrysene
- 12 Benzo[b]fluoranthene
- 13 Benzo[k]fluoranthene
- 14 Benzo[a]pyrene
- 15 Dibenzo[a,h]anthracene
- 16 Indeno[1,2,3-cd]pyrene
- 17 Benzo[g,h,i]perylene

