

Soft Drink Additives on Kinetex 2.6µm C18 150 x 4.6mm ID

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Gradient

Eluent A: 0.1% phosphoric acid in water


Eluent B: 0.1% phosphoric acid in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1.5	95	5
	3	6	60	40
	4	9.75	60	40
	5	11.25	5	95
	6	11.26	95	5
	7	15.75	95	5

Flow Rate: 1.2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 215 nm (ambient)

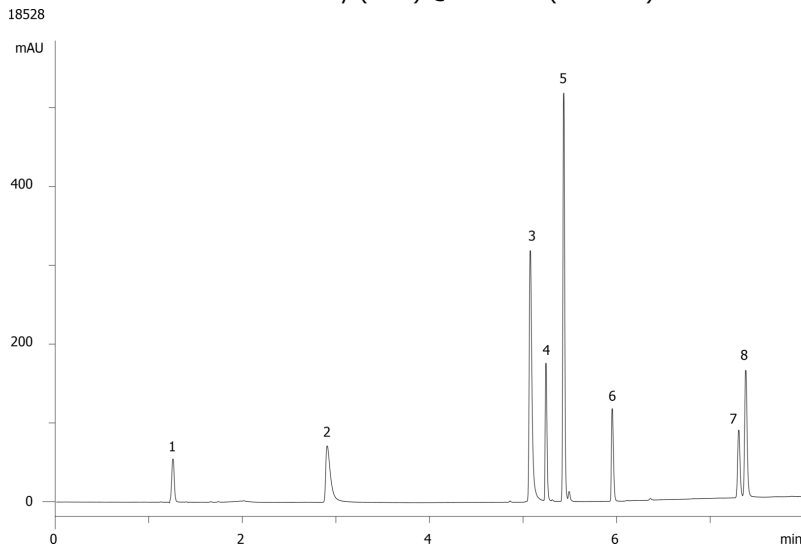

ANALYTES:

- 1** Ascorbic acid
Retention Time: 1.258 min
- 2** Acesulfame-K
Retention Time: 2.907 min
- 3** Saccharin
Retention Time: 5.077 min
- 4** Quinine
Retention Time: 5.245 min
- 5** Caffeine
Retention Time: 5.435 min
- 6** Aspartame
Retention Time: 5.953 min
- 7** Sorbic acid
Retention Time: 7.305 min
- 8** Benzoic acid
Retention Time: 7.38 min

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com



Phenomenex products are available worldwide.

www.phenomenex.com.cn

info@phenomenex.com