

QC Test Chromatogram (Agilent 1100) - Kinetex 2.6µm C18 50 x 4.6mm ID

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4462-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50


Flow Rate: 2.35 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

18565

4

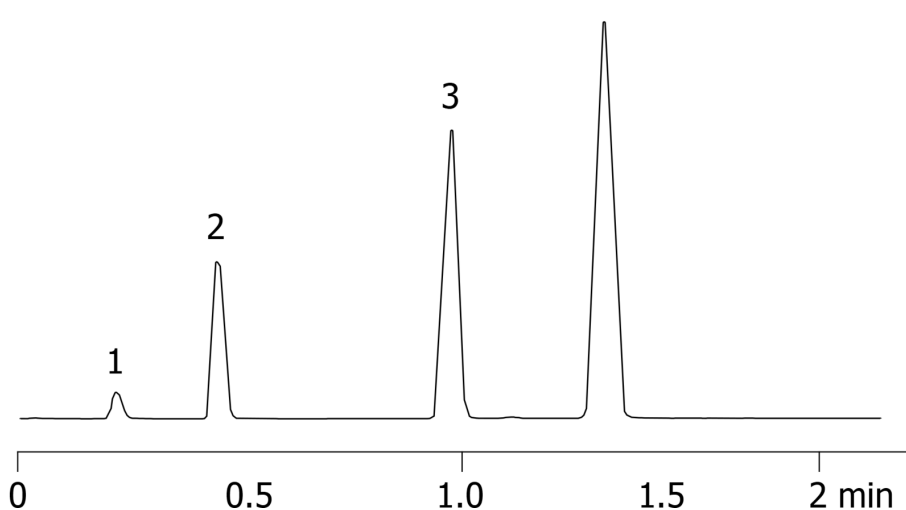

ANALYTES:

- 1** Uracil
Retention Time: 0.22 min
- 2** Acetophenone
Retention Time: 0.45 min
- 3** Toluene
Retention Time: 0.98 min
- 4** Naphthalene
Retention Time: 1.32 min



Products used in this application:

