

Epinephrine and Norepinephrine on Kinetex 2.6µm HILIC 50 x 2.1mm ID

Column: Kinetex® 2.6 µm HILIC 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4461-AN

Elution Type: Isocratic

Eluent A: Acetonitrile:

Eluent B: 100 mM Ammonium formate pH 3.2

Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	92	8

Profile:

Flow Rate: 0.4 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

18605



Products used in this application:

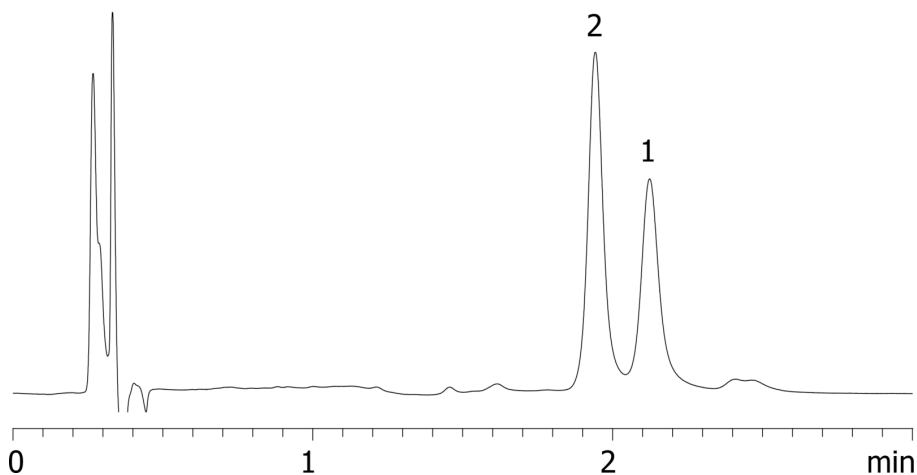

ANALYTES:

- 1 Norpinephrine
Retention Time: 2.123 min
- 2 Epinephrine
Retention Time: 1.942 min

