

QC Test Chromatogram (Agilent 1200SL) - Luna 3µm C18(2) 150 x 4.6mm ID

Column: Luna® 3 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4251-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30

Flow Rate: 1.8 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Detector Info: Agilent 1200SL. BP = 3625 psi; N = 166,502 p/m


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

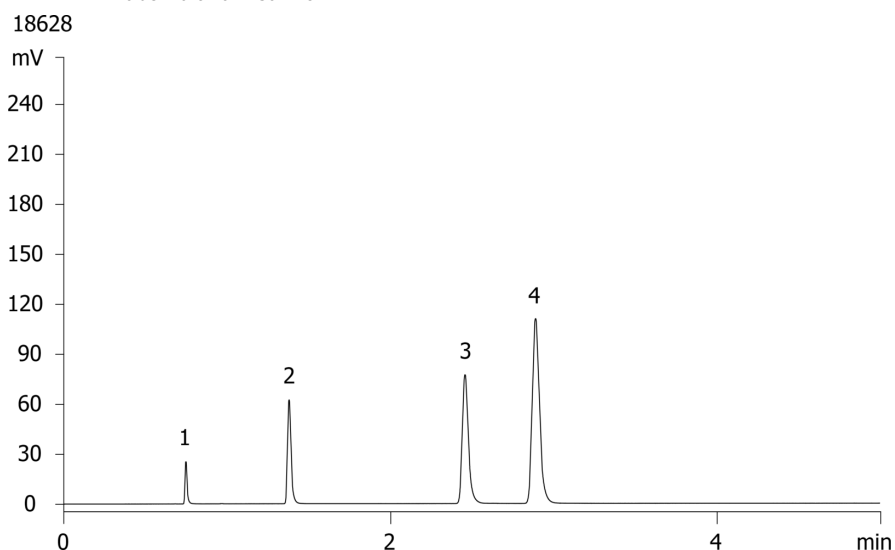

ANALYTES:

- 1** Uracil
Retention Time: 0.75 min
- 2** Acetophenone
Retention Time: 1.38 min
- 3** Toluene
Retention Time: 2.46 min
- 4** Naphthalene
Retention Time: 2.89 min

