

Extraction of PCP using Strata X-C on Kinetex 2.6µm C18 50 x 2.1mm ID with Acquity UPLC MSMS

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% Formic acid in H2O

Eluent B: 0.1% Formic acid in ACN


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	0.25	97	3
	3	2.75	5	95
	4	3	5	95
	5	3.01	97	3
	6	4.5	97	3

Flow Rate: 500 µL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (110 °C)

18652

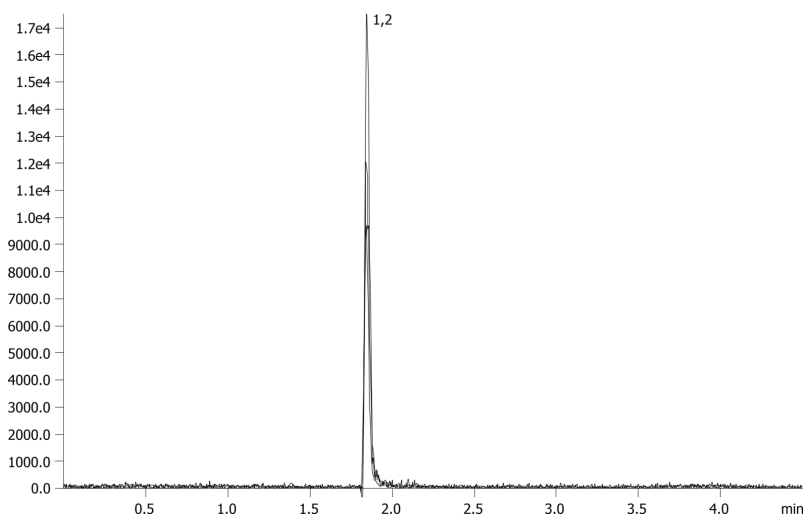

ANALYTES:

1 PCP

Retention Time: 1.85 min

2 PCP-d5

Retention Time: 1.84 min



Products used in this application:

