

Steroids on Kinetex 2.6µm C18 150 x 4.6mm ID

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: CH₃CN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	12	35	65
	3	12.1	75	25
	4	16	75	25


Flow Rate: 1.2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 230 nm (ambient)



Products used in this application:

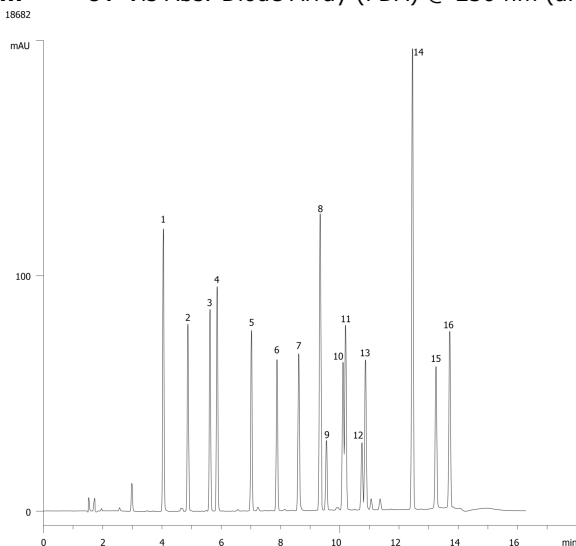

Steroids on Kinetex 2.6µm C18 150 x 4.6mm ID

ANALYTES:

- 1** Triamcinolone
- 2** Estriol
- 3** Prednisone
- 4** Cortisone
- 5** Betamethasone
- 6** Corticosterone
- 7** Triamcinalone acetate
- 8** 11-alpha-Hydroxyprogesterone
- 9** Cortisone acetate
- 10** Estradiol
- 11** 21-Hydroxyprogesterone
- 12** 11-Ketoprogesterone
- 13** Estrone
- 14** 17-Hydroxyprogesterone
- 15** Deoxycorticosterone
- 16** Betamethasone Valerate

