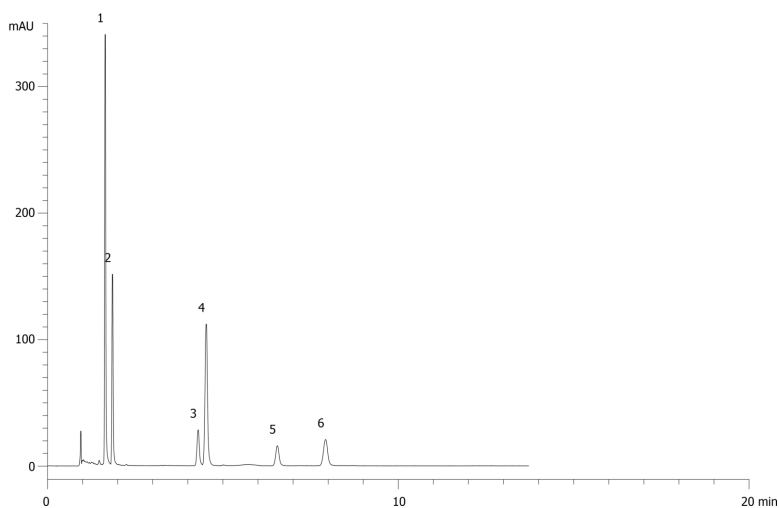


Fat Soluble Vitamins on Kinetex 2.6µm C18 150 x 4.6mm ID**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4462-E0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Eluent B:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5

Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

18696

**Products used in this application:****ANALYTES:**

- 1** Vitamin A
Retention Time: 1.64 min
- 2** Vitamin A acetate
Retention Time: 1.85 min
- 3** Vitamin E
Retention Time: 4.29 min
- 4** Vitamin D3
Retention Time: 4.52 min
- 5** Vitamin E acetate
Retention Time: 6.55 min
- 6** Vitamin K1
Retention Time: 7.92 min

