

Prodigy Test Mix on Kinetex 2.6um C18 50 x 3.0 mm ID on HP1100

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 3 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4462-Y0

Elution Type: Isocratic

Eluent A: Acetonitrile


Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

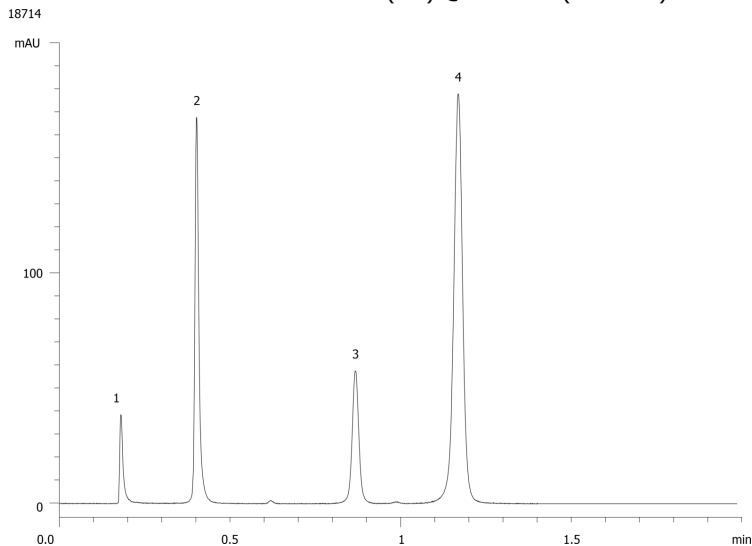

ANALYTES:

- 1 Uracil
- 2 Acetophenone
- 3 Toluene
- 4 Naphthalene



Products used in this application:

