

FMOC-DL-Tryptophan-(boc)-OH on Lux 5 µm Cellulose-2 250 x 4.6mm in NP**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1% TFA in Hexane**Eluent B:** 0.1% TFA in IPA

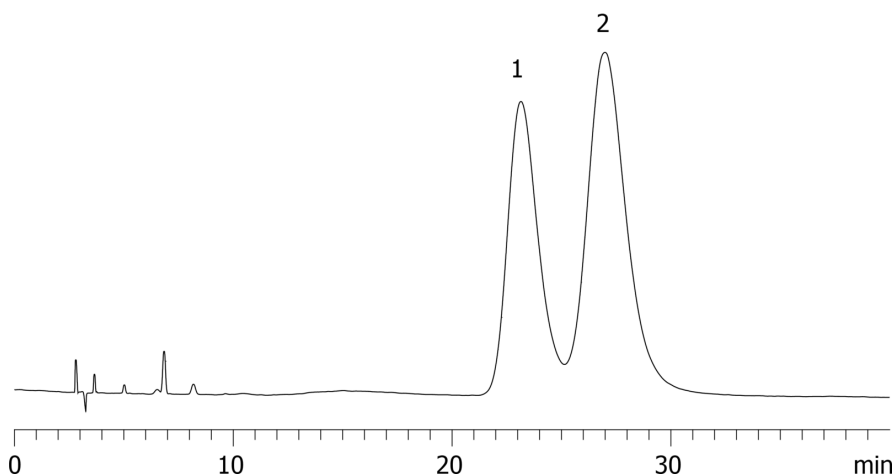
Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	30	90	10

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

18733

**Products used in this application:****ANALYTES:****1** FMOC-Trp(Boc)-OH

Retention Time: 23.1 min

2 FMOC-Trp(Boc)-OH

Retention Time: 26.9 min

