

HPLC Application

ID No.: 18757

Kinetex 2.6µm C18 Batch 5569-68

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4462-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	65	35

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

18757



Products used in this application:

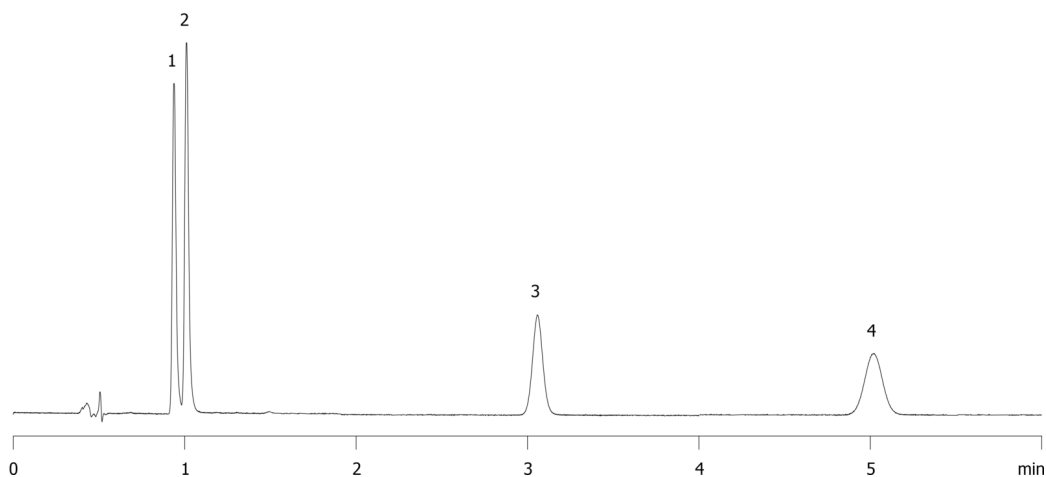

ANALYTES:

- 1 Uracil
- 2 Hydroxycortisone
- 3 Cortisone
- 4 Cortisone acetate
- 5 17-Hydroxycortisone

