

Pharma Screen in 3.77min using Kinetex 2.6u C18 100x4.6mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4462-E0

Elution Type: Gradient

Eluent A: Water with 0.1% H3PO4

Eluent B: Acetonitrile with 0.1% H3PO4

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	5	95

Flow Rate: 2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 215 nm (ambient)

18758



Products used in this application:

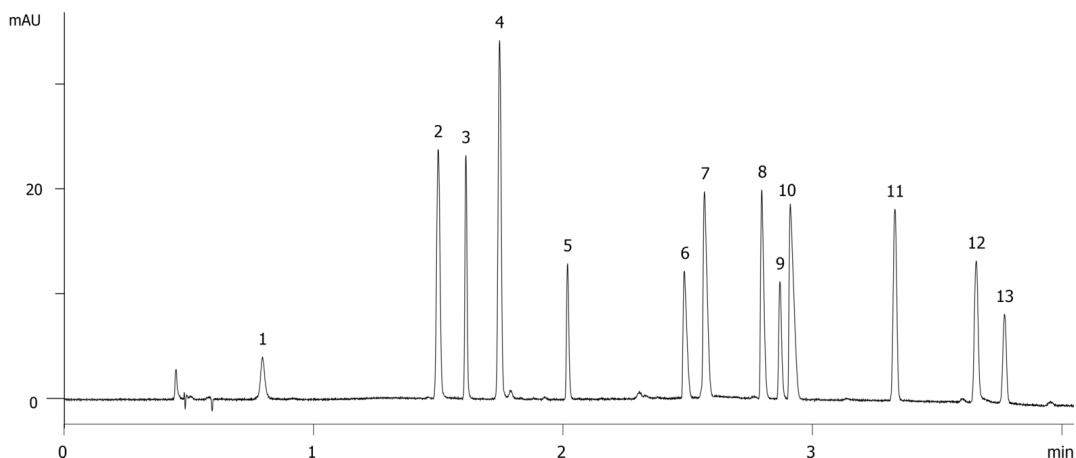

Pharma Screen in 3.77min using Kinetex 2.6u C18 100x4.6mm

ANALYTES:

- 1** Procainamide
- 2** Acetaminophen
- 3** Folic acid
- 4** Sulfathiazole
- 5** Acebutolol
- 6** Dextromethorphan
- 7** Diphenhydramine
- 8** Propafenone
- 9** Amitriptyline
- 10** Fluoxetine
- 11** Naproxen
- 12** Diflunisal
- 13** Indomethacin

