


EPA Method 8081- Organochlorine Pesticides on ZB-35HT**Column:** Zebron[™] ZB-35HT, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G025-11**Oven Profile:** 160 °C for 0.5 min to 320 °C @ 25 °C/min for 1.5 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 5:1 1 μ L @ 350°C**Detection:** Electron Capture (ECD) (350°C)**Analyst Note:** Analytes are 50ppb in isooctane**Products used in this application:**

EPA Method 8081- Organochlorine Pesticides on ZB-35HT

ANALYTES:

- 1 alpha-BHC
- 2 gamma-BHC
- 3 b-BHC
- 4 Heptachlor
- 5 delta-BHC
- 6 Aldrin
- 7 Heptachlor epoxide
- 8 gamma-Chlordane
- 9 alpha-Chlordane
- 10 Endosulfan I
- 11 DDE
- 12 Dieldrin
- 13 Endrin
- 14 DDD
- 15 Endosulfan II
- 16 DDT
- 17 Endrin aldehyde
- 18 Endosulfan sulfate
- 19 Methoxychlor
- 20 Endrin ketone

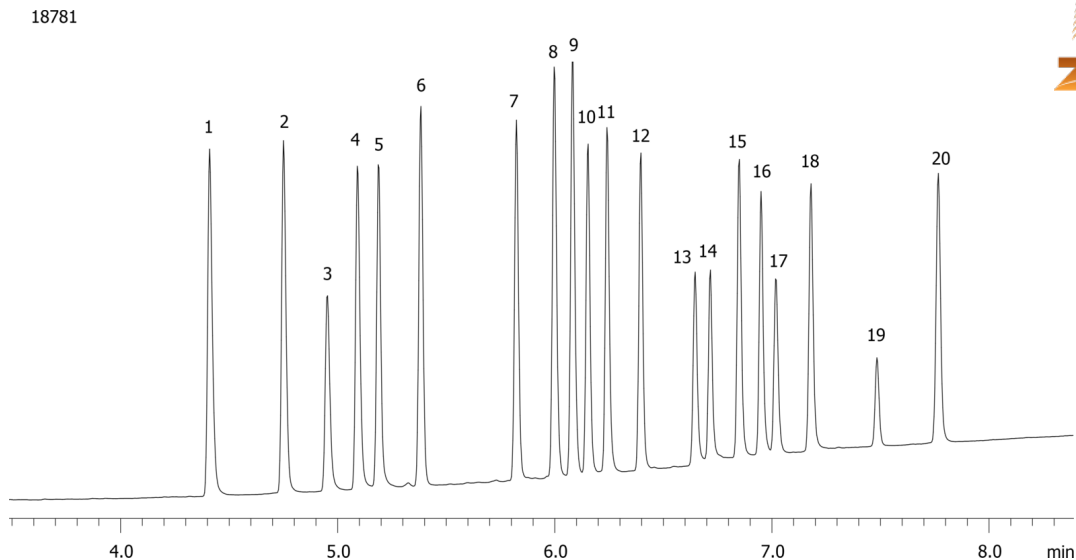