


EPA Method 8081 - Chlorinated Pesticides on ZB-XLB-HT**Column:** Zebron™ ZB-XLB-HT, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G024-11**Oven Profile:** 110 °C to 320 °C @ 15 °C/min for 2 min**Carrier Gas:** Constant Flow Helium, 0.65 mL/min**Injection:** Split 15:1 1 µL @ 250°C**Detection:** Electron Capture (ECD) (350°C)**Analyst Note:** Analytes are 20ppb in isooctane.**Products used in this application:**

EPA Method 8081 - Chlorinated Pesticides on ZB-XLB-HT

ANALYTES:

- 1 alpha-BHC
- 2 gamma-BHC
- 3 b-BHC
- 4 delta-BHC
- 5 Heptachlor
- 6 Aldrin
- 7 Heptachlor epoxide
- 8 gamma-Chlordane
- 9 alpha-Chlordane
- 10 Endosulfan I
- 11 4,4'-DDE
- 12 Dieldrin
- 13 Endrin
- 14 4,4'-DDD
- 15 Endosulfan II
- 16 Endrin aldehyde
- 17 4,4'-DDT
- 18 Endosulfan sulfate
- 19 Methoxychlor
- 20 Endrin ketone

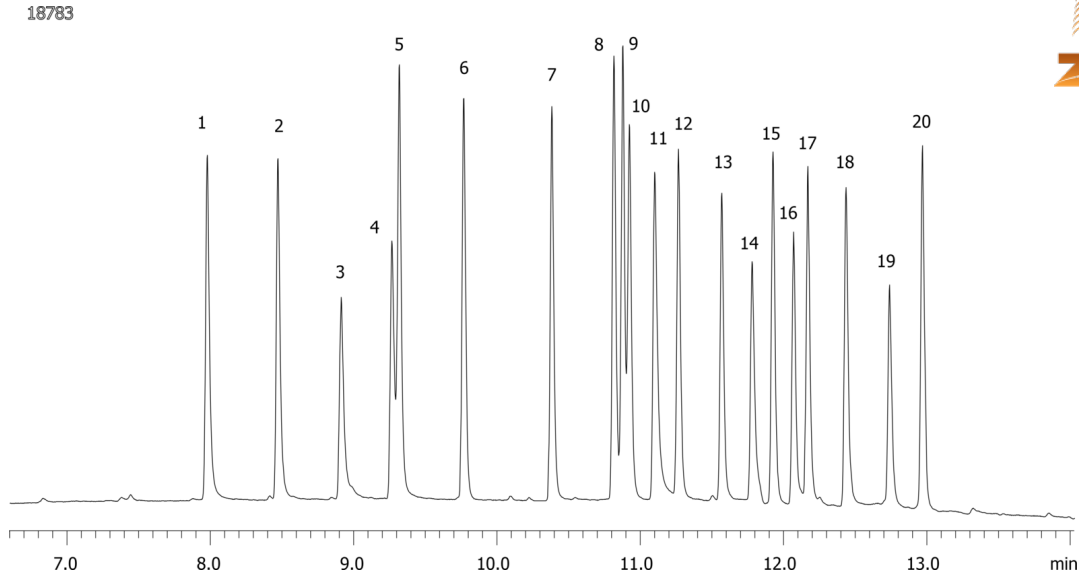