

SecurityGuard Comparisons on Gemini 5u C18 150 x 4.6mm With Cartridge

Column: Gemini® 5 µm C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35

Flow Rate: 1

Col. Temp.: 25 °C


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: Analysis performed with SecurityGuard Cartridge PN: AJ0-7597 (cf. App IDs: 18795 and 18797) using Prodigy Test Mix (Cat #: 0-3045)

18796



Products used in this application:

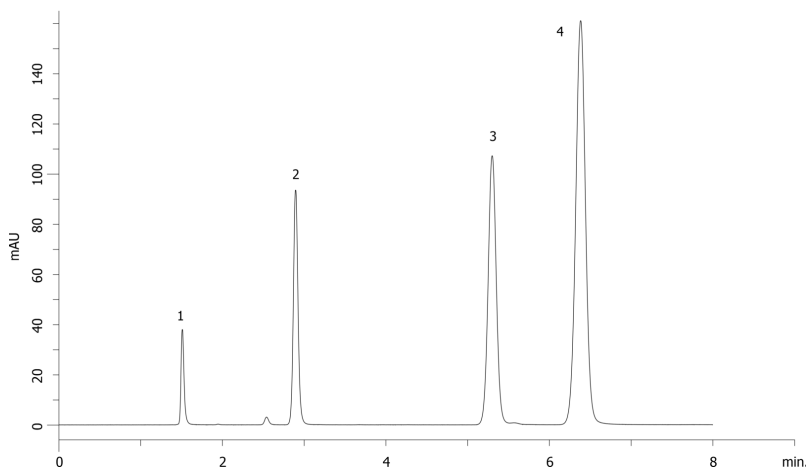

ANALYTES:

1 Uracil

Retention Time: 1.51 min

2 Acetophenone

Retention Time: 2.895 min

3 Toluene

Retention Time: 5.3 min

4 Naphthalene

Retention Time: 6.382 min

