

SecurityGuard Comparisons on Luna 5u C18(2) 150 x 4.6mm Without Cartridge

Column: Luna® 5 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	65	35	

Flow Rate: 1 mL/min

Col. Temp.: 25 °C


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: Analysis performed without SecurityGuard Cartridge (cf. App IDs: 18808 and 18809) using Prodigy Test Mix (PN: AL0-3045)

18807



Products used in this application:

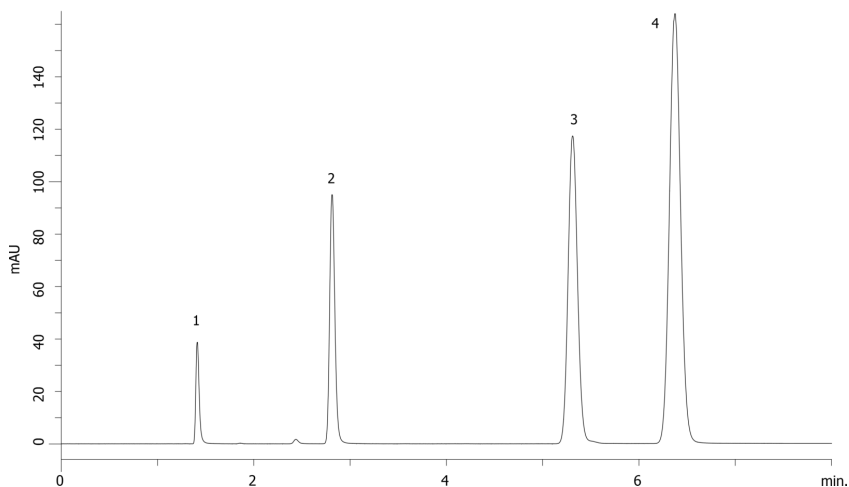

ANALYTES:

- 1** Uracil
Retention Time: 1.413 min
- 2** Acetophenone
Retention Time: 2.813 min
- 3** Toluene
Retention Time: 5.311 min
- 4** Naphthalene
Retention Time: 6.372 min

