

Atenolol and Related Substances CRS using Luna 5µm C18(2) 125 x 4.0mm

Column: Luna® 5 µm C18(2) 100 Å, LC Column 125 x 4 mm, Ea

Dimensions: 125 x 4 mm ID

Order No: 00E-4252-D0

Elution Type: Isocratic

Eluent A: EP Mobile Phase

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.6 mL/min


Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 226 nm (ambient)

Analyst Note: 12.5 mM KH₂PO₄ in Water pH 3.0 + 2.0g Sodium Octanesulfonate + 0.8g Tetrabutyl Ammonium Hydrogen Sulfate/Methanol/THF v/v/v 80:18:2
18825



Products used in this application:

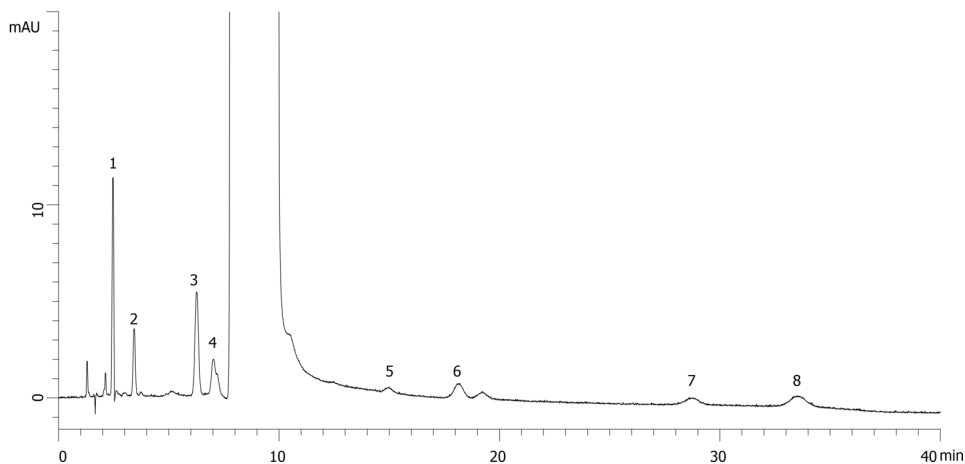

ANALYTES:

- 1 Impurity B
- 2 Impurity A
- 3 Impurity J
- 4 Impurity I
- 5 impurities D and E
- 6 Impurity F
- 7 Impurity G
- 8 Impurity H

