

Extraction of Expanded Amphetamines Using Strata X-C on Kinetex 2.6µm C18 50 x 2.1mm ID by LCMS

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1%FA-5mM ammonium acetate in H2O

Eluent B: 0.1%FA-5mM ammonium acetate in Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	3	60	40
	3	3.02	5	95
	4	3.3	5	95
	5	3.31	90	10
	6	5.1	90	10

Flow Rate: 450 µL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)

Detector Info: ESI conditions

Positive Ionization


NEB:	12.00
CUR:	11.00
IS:	4500.00
TEM:	450.00
CAD:	11.00
EP:	10.00



Products used in this application:



18847

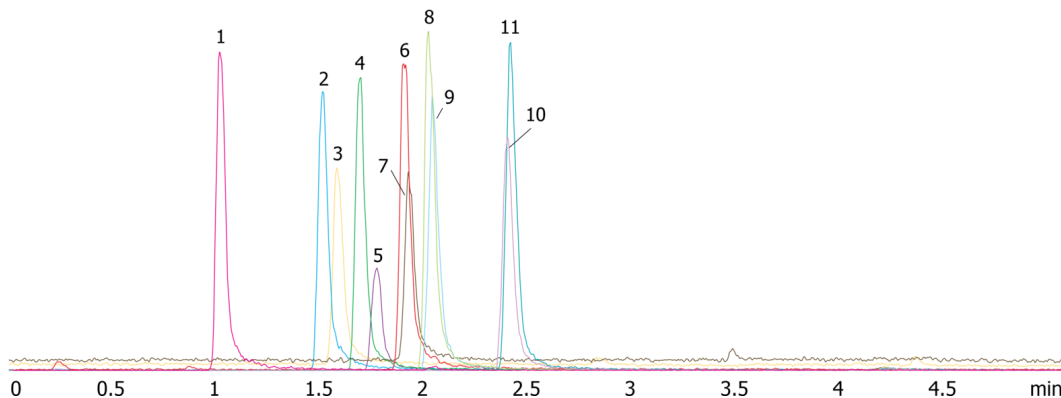

Extraction of Expanded Amphetamines Using Strata X-C on Kinetex 2.6µm C18 50 x 2.1mm ID by LCMS

ANALYTES:

- 1** HMMA
- 2** Amphetamine-d11
- 3** Amphetamine
- 4** Methamphetamine
- 5** Methamphetamine-d14
- 6** MDA-d5
- 7** MDA
- 8** MDMA-d5
- 9** MDMA
- 10** MDEA-d5
- 11** MDEA

