


Mirtazapine on Lux 5 µm Cellulose-2 by Chiral SFC

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4457-E0
Elution Type: Isocratic
Eluent A: Methanol with 0.1% DEA
Eluent B: CO₂
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	20	80

Profile:
Flow Rate: 2.5 mL/min
Col. Temp.: 35 °C
Detection: UV-Vis Abs.-Diode Array (PDA) @ 220 nm (ambient)
Analyst Note: Application using Jasco SFC instrument courtesy of Cephalon, Inc.

18856

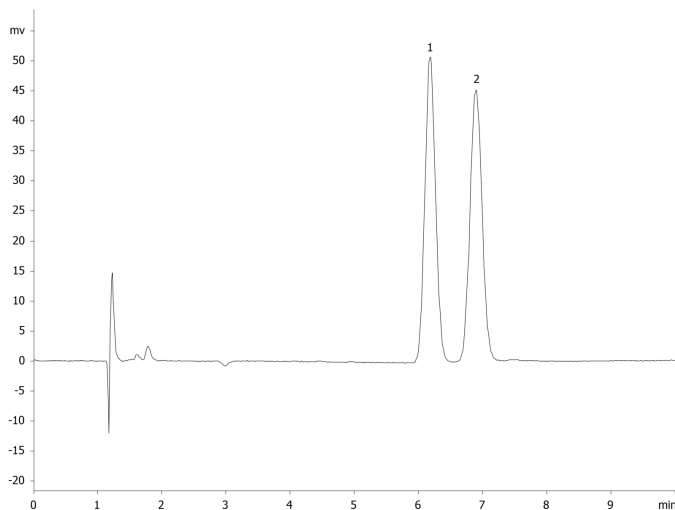

ANALYTES:

- 1 Mirtazapine
- 2 Mirtazapine



Products used in this application:

