

Bopindolol on Lux 5 µm Cellulose-1 by Chiral SFC

Column: Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4459-E0

Elution Type: Isocratic

Eluent A: Ethanol with 0.1% DEA

Eluent B: CO₂


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	0	30	70

Flow Rate: 2.5 mL/min

Col. Temp.: 35 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 220 nm (ambient)

Analyst Note: Application using Jasco SFC instrument courtesy of Cephalon, Inc.

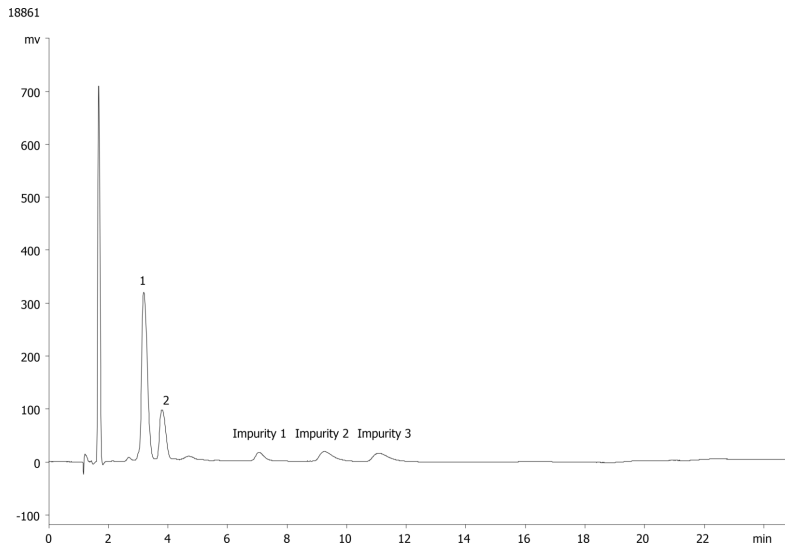

ANALYTES:

- 1 Bopindolol
- 2 Bopindolol



Products used in this application:

