

Mycotoxins on Kinetex 2.6µm C18 50 x 2.1mm ID by LCMS in ESI Positive Mode (1)

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% Formic acid

Eluent B: 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	3	5	95
	3	4	5	95
	4	4.01	70	30
	5	6	70	30

Flow Rate: 400 µL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)

Detector Info: ESI conditions


Positive Ionization

NEB:	12.00
CUR:	11.00
IS:	4500.00
TEM:	450.00
CAD:	11.00
EP:	10.00

App ID 18878



Products used in this application:

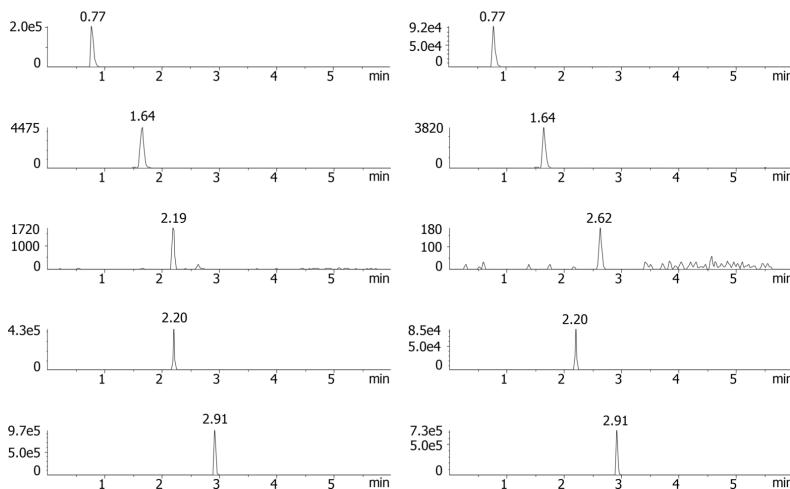

Mycotoxins on Kinetex 2.6µm C18 50 x 2.1mm ID by LCMS in ESI Positive Mode (1)

ANALYTES:

- 1** Fumonisin B1
- 2** HT-2
- 3** T2-Toxin
- 4** Fumonisin B2
- 5** Ochratoxin

