

PEGylated IgG (N-Terminal PEG 40KDa) on BioSep3000 (2)

Column: BioSep™ 5 µm SEC-s3000 290 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2146-K0

Elution Type: Isocratic

Eluent A: 100mM Sodium Phosphate pH 6.8□□

Gradient **Step No.** **Time (min)** **Pct A**

Profile: **1** 0 100

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, GFC 3000 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4488


- Holder Part No.: KJ0-4282

18939

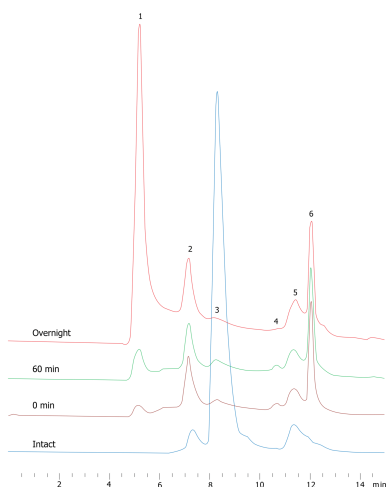

ANALYTES:

- 1 High MW PEG/IgG complex
- 2 IgG dimer / IgG / 1 PEG complex
- 3 IgG monomer
- 4 impurity
- 5 impurity
- 6 impurity

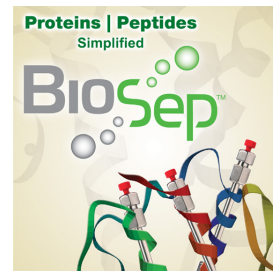

Products used in this application:

