

PEGylated alpha-Chymotrypsinogen A (N-Terminal PEG 20KDa) on BioSep3000 (2)

Column: BioSep™ 5 µm SEC-s3000 290 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2146-K0

Elution Type: Isocratic

Eluent A: 100mM Sodium Phosphate pH 6.8□□

Gradient **Step No.** **Time (min)** **Pct A**

Profile: **1** 0 100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, GFC 3000 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4488

- Holder Part No.: KJ0-4282

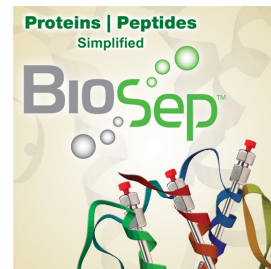

Products used in this application:

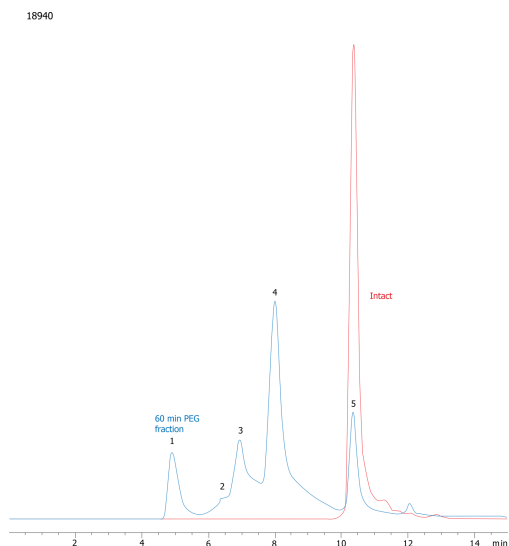

ANALYTES:

- 1 PEG/Protein
- 2 3 PEG / Protein complex
- 3 2 PEG / Protein complex
- 4 PEGylated Chymotrypsinogen A
- 5 Chymotrypsinogen A

