

Indapamide on Lux 5 µm Cellulose-1 by Chiral SFC


Column: Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4459-E0
Elution Type: Isocratic
Eluent A: Methanol with 0.1% DEA
Eluent B: Carbon Dioxide
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	20	80

Profile:
Flow Rate: 1 mL/min
Col. Temp.: 35 °C
Detection: UV-Vis Abs.-Diode Array (PDA) @ 220 nm (ambient)
Analyst Note: Application using Jasco SFC instrument courtesy of Cephalon, Inc.
 18951



Products used in this application:

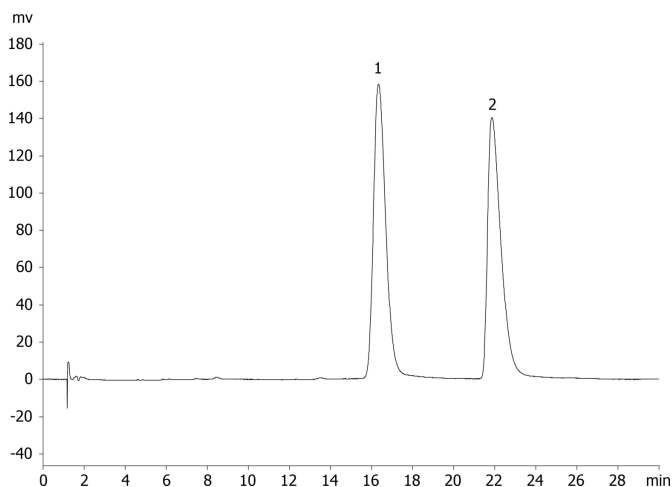

ANALYTES:

- 1 Indapamide
- 2 Indapamide

