

Glibenclamide Related Substance HPLC Analysis (HP-EUR) on Kinetex 2.6u C18 50 x 4.6mm ID

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4462-E0

Elution Type: Gradient

Eluent A: 20:50:930 101.8g/L Triethylamine solution/ACN/H2O


Eluent B: 20:65:915 mobile phase A/H2O/ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	45	55
	2	15	45	55
	3	30	5	95
	4	40	5	95
	5	41	45	55
	6	55	45	55

Flow Rate: 0.8 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)

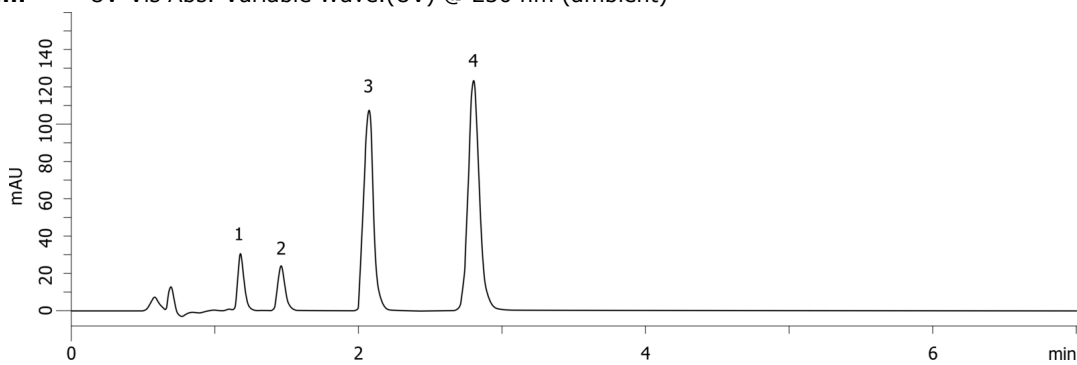

Products used in this application:



ANALYTES:

1 Glibenclamide impurity A

Retention Time: 1.18 min

2 Glibenclamide impurity B

Retention Time: 1.46 min

3 Glibenclamide

Retention Time: 2.068 min

4 Glizlazide

Retention Time: 2.799 min

