

SecurityGuard Comparisons on Lux 250 x 4.6mm Without Cartridge

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Acetonitrile

Eluent B: 20mM NH₄HCO₃


Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	60	40

Flow Rate: 0.7 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: Analysis performed without SecurityGuard Cartridge (cf. App IDs: 18823 and 18824) using Prodigy Test Mix (PN: AL0-3045) 19070

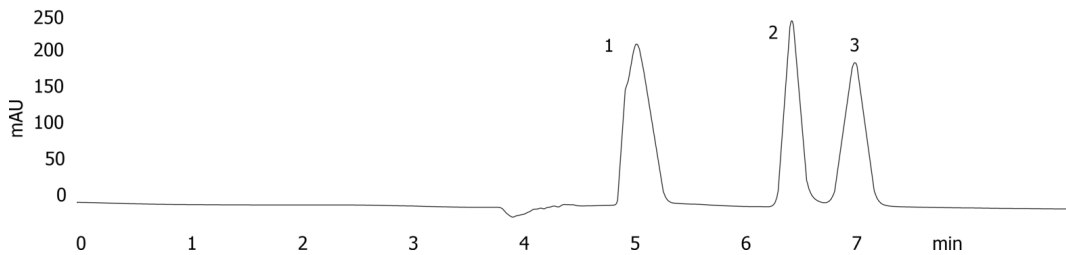

ANALYTES:

- 1 impurity
- 2 Chlorpheniramine
- 3 Chlorpheniramine

