

Ginsenoside standards using Kinetex 2.6u C18 100x4.6mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4462-E0

Elution Type: Gradient

Eluent A: Water with 0.1% H₃PO₄

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	78	22
	2	10	78	22
	3	20	50	50

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 200 nm (ambient)

19074



Products used in this application:

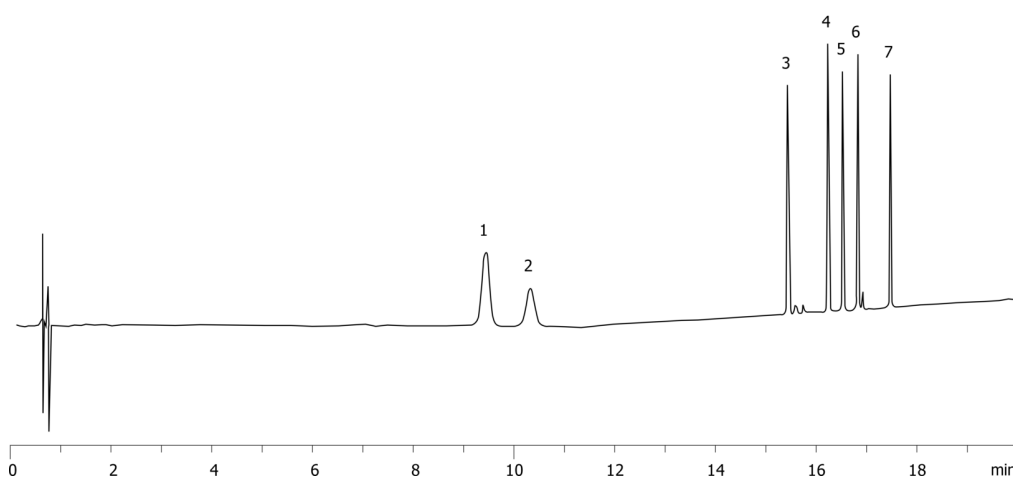

ANALYTES:

- 1** Ginsenoside Rg1
- 2** Ginsenoside Re
- 3** Ginsenoside Rf
- 4** Ginsenoside Rb1
- 5** Ginsenoside Rc
- 6** Ginsenoside Rb2
- 7** Ginsenoside Rg

