

Catechins on Kinetex 2.6 µm C18 100 x 4.6 mm at 1.5 ml/min

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4462-E0

Elution Type: Gradient

Eluent A: 0.1% Formica acid

Eluent B: 0.1% Formic acid in CH₃CN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	0.53	85	15
	3	3.73	65	35
	4	3.74	85	15
	5	5.87	85	15

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 230 nm (ambient)

19075



Products used in this application:

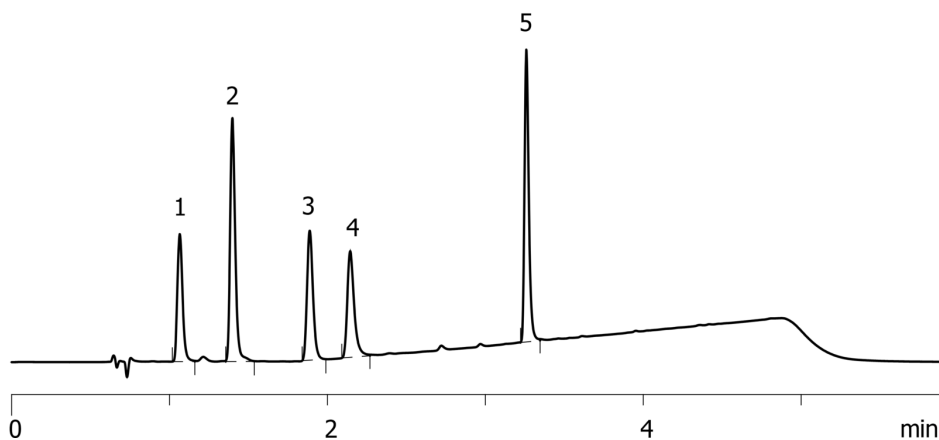

ANALYTES:

- 1** Epigallocatechin
Retention Time: 1.065 min
- 2** Catechin
Retention Time: 1.398 min
- 3** Epicatechin
Retention Time: 1.888 min
- 4** Epigallocatechin gallate
Retention Time: 2.144 min
- 5** Epicatechin gallate
Retention Time: 3.26 min

