

Compare BioSep-SEC-S2000(2) & Sepax SRT SEC-150 4 proteins

Column: BioSep™ 5 µm SEC-s2000 145 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2145-K0

Elution Type: Isocratic

Eluent A: 10mM Tris, 150mM NaCl, pH 7.4

Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

Flow Rate: 0.6 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, GFC-2000 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4487


- Holder Part No.: KJ0-4282

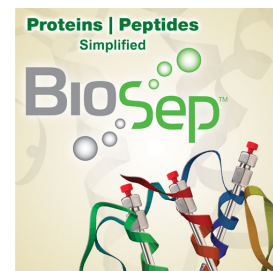

Products used in this application:



19091

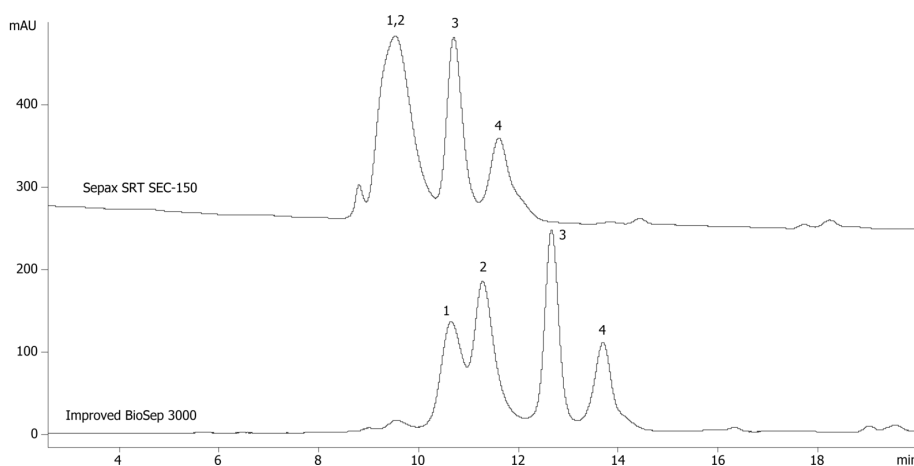

ANALYTES:

- 1 IgA
- 2 beta-Amylase
- 3 BSA
- 4 Ovalbumin

