

Comparison of Cellulose 1, 2 and Amylose 2 Selectivity with IPA**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Isopropyl Alcohol with 0.1% DEA**Eluent B:** Carbon Dioxide

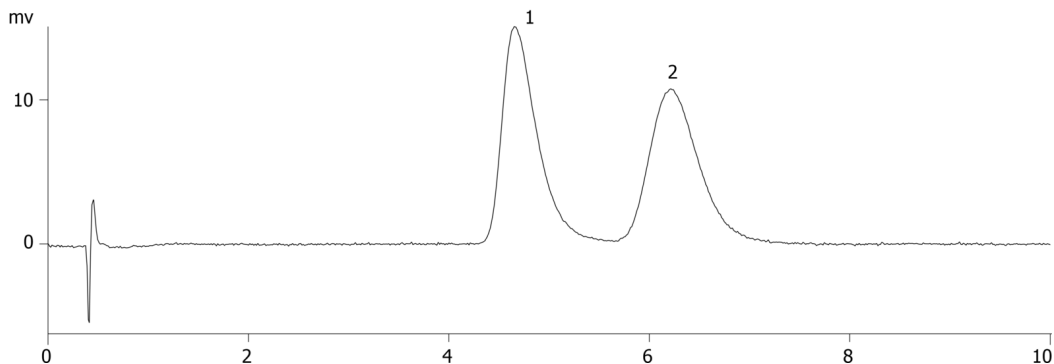
Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	20	80

Flow Rate: 2.5 mL/min**Col. Temp.:** 35 °C**Detection:** UV-Vis Abs.-Fixed Wave. (UV) @ 254 nm (35 °C)**Detector Info:** PDA**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Silica 4 x 3.0mm, 10/Pk Part No.: AJ0-4348

- Holder Part No.: KJ0-4282

19110

**ANALYTES:****1** Propafenone