

Lysozyme analysis using Kinetex 2.6u C18 50x4.6mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4462-E0

Elution Type: Gradient

Eluent A: Water with 0.1% TFA


Eluent B: Acetonitrile with 0.1% TFA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	5	5	95

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

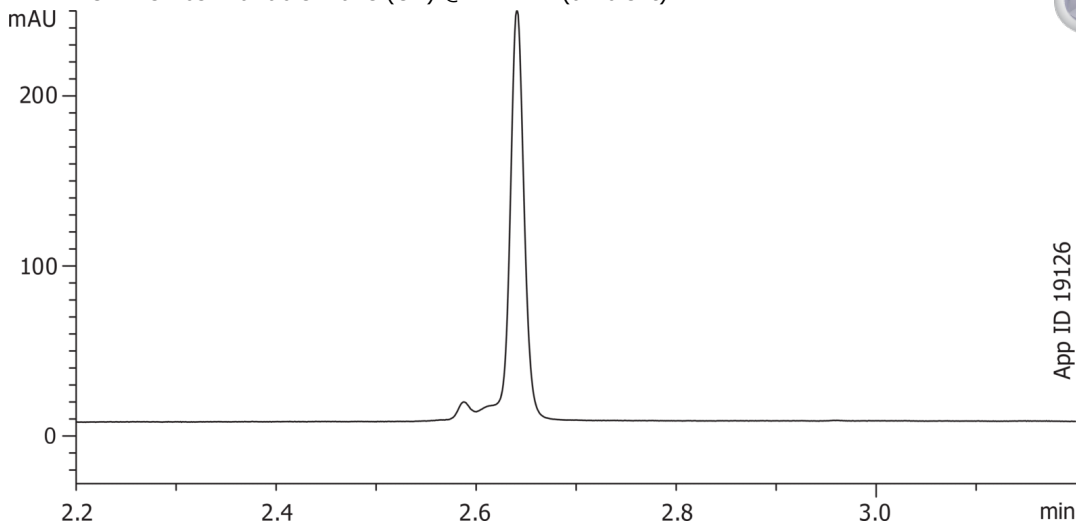

Products used in this application:



ANALYTES:

1 Lysozyme

