

HPLC Application

ID No.: 19134

Hu IgG2 k Lys-C & Asp-N Digest on Kinetex C18 2.6u

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Gradient

Eluent A: 0.1% TFA and 2% Acetonitrile in Water

Eluent B: 0.085% TFA in Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	2	98	2
	2	35	52	48
	3	40	30	70


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (25 °C)
19134



Products used in this application:

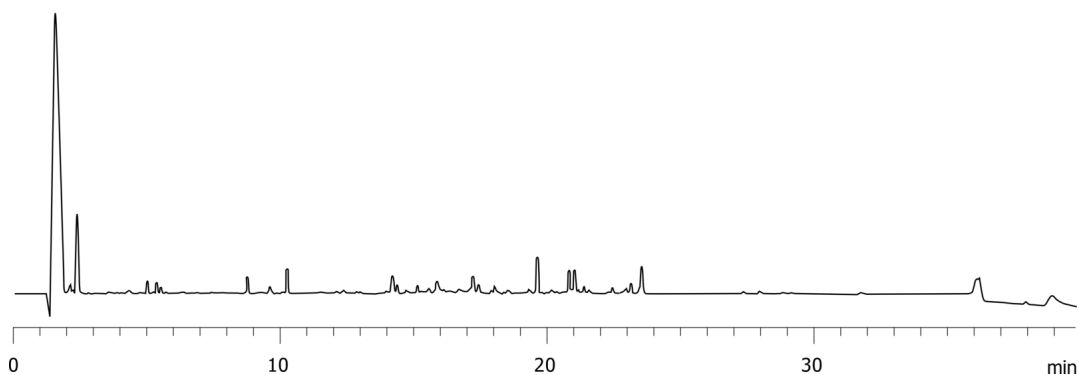

ANALYTES:

- 1 Human IgG2 k Lys-C & Asp-N