

QC Test Chromatogram (Shimadzu Prominence UFLCXR) - Kinetex 2.6µm C18 50 x 3mm ID


Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 3 mm, Ea
Dimensions: 50 x 3 mm ID
Order No: 00B-4462-Y0
Elution Type: Isocratic
Eluent A: Acetonitrile
Eluent B: Water
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	50	50

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (25 °C)
 19145



Products used in this application:

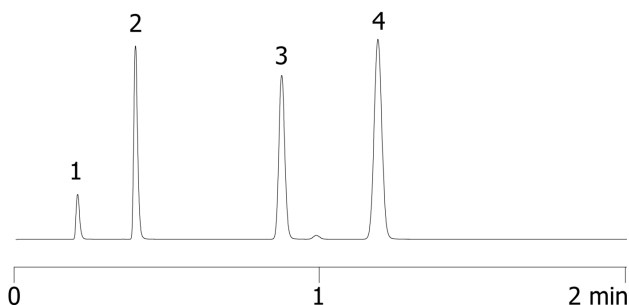

ANALYTES:

- 1** Uracil
Retention Time: 0.2 min
- 2** Acetophenone
Retention Time: 0.435 min
- 3** Toluene
Retention Time: 0.951 min
- 4** Naphthalene
Retention Time: 1.296 min

