

Poly dT12-18 Oligonucleotides on Luna HILIC

Column: Luna® 3 µm HILIC 200 Å, LC Column 50 x 2 mm, Ea

Dimensions: 50 x 2 mm ID

Order No: 00B-4449-B0

Elution Type: Gradient

Eluent A: AcCN

Eluent B: 10 mM NH₄OAc

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	10	50	50

Flow Rate: 0.3 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, HILIC 4 x 2.0mm ID, 10/Pk Part No.: AJ0-8328


- Holder Part No.: KJ0-4282



Products used in this application:



19176

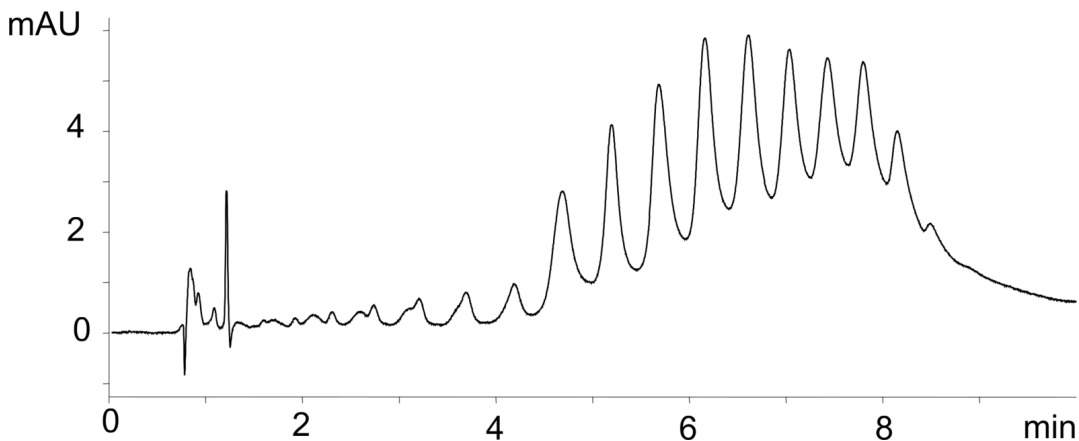

ANALYTES:

1 poly dT 12-18 mixture

