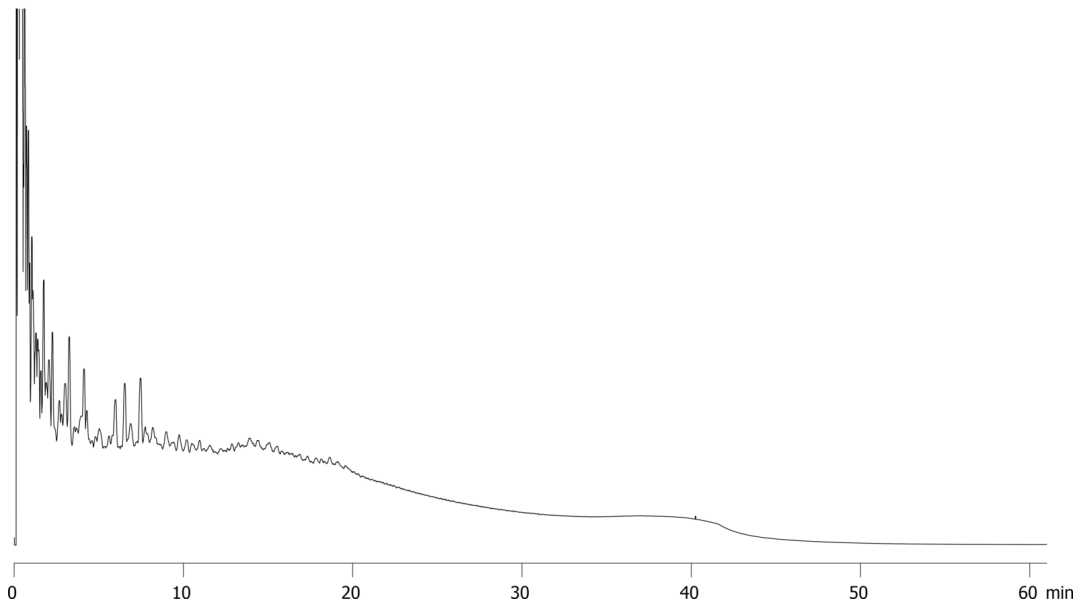


Yeates Crude Oil on ZB-1XT 5x0.53x0.09

Column: Zebron™ ZB-1XT SimDist Metal, GC Cap. Column 5 m x 0.53 mm x 0.09 µm, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 5 meters x 0.53 mm x 0.09 µm
Order No: 7AK-G026-55
Oven Profile: 35 °C to 450 °C @ 10 °C/min for 20 min

Carrier Gas: Constant Flow Helium, 15 mL/min
Injection: On-Column :1 1 µL @ 38°C
Detection: Flame Ionization (FID) (450°C)
Analyst Note: Sample is 1% in CS₂. Chromatogram is baseline subtracted.

19256

**ANALYTES:**

- 1 Yeates crude oil

**Products used in this application:**