


**ASTM D7169 (Hydrocarbons C7-C100+) on ZB-1XT 5x0.53x0.15****Column:** Zebtron<sup>™</sup> ZB-1XT SimDist Metal, GC Cap. Column 5 m x 0.53 mm x 0.15  $\mu$ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 5 meters x 0.53 mm x 0.15  $\mu$ m**Order No:** 7AK-G026-05**Oven Profile:** 30 °C to 450 °C @ 10 °C/min for 10 min**Carrier Gas:** Constant Flow Helium, 7 mL/min**Injection:** On-Column :1 1  $\mu$ L @ 33°C**Detection:** Flame Ionization (FID) (450°C)**Analyst Note:** Sample in CS<sub>2</sub>. Chromatogram is baseline subtracted.**Products used in this application:**

19264

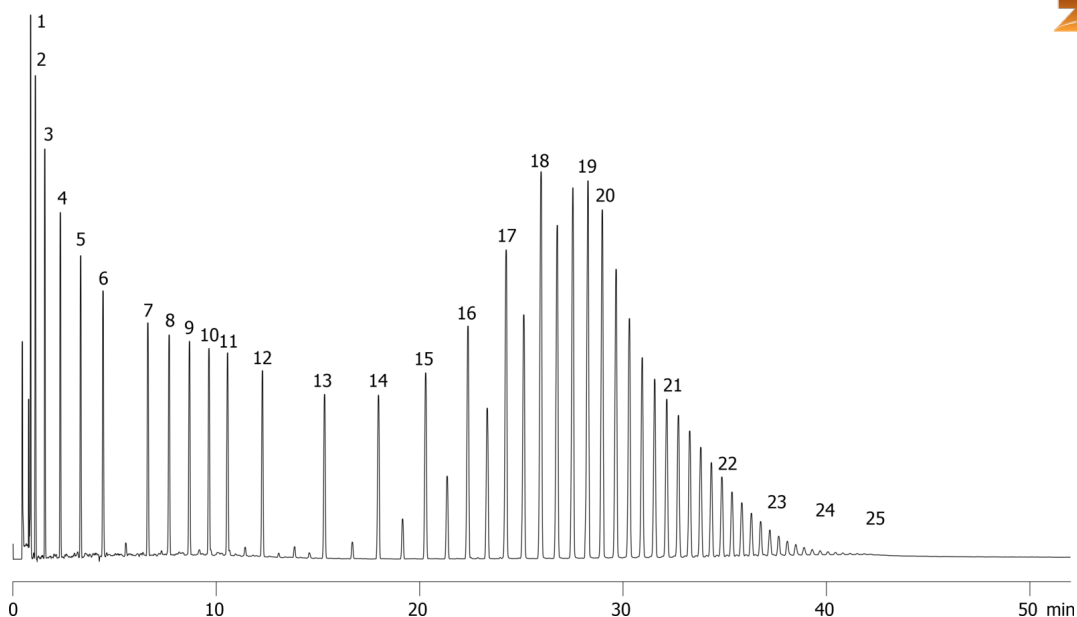

## ASTM D7169 (Hydrocarbons C7-C100+) on ZB-1XT 5x0.53x0.15

### ANALYTES:

- 1 Heptane (C7)
- 2 Octane (C8)
- 3 Nonane (C9)
- 4 Decane (C10)
- 5 Undecane (C11)
- 6 Dodecane (C12)
- 7 Tetradecane (C14)
- 8 Pentadecane (C15)
- 9 Hexadecane (C16)
- 10 Heptadecane (C17)
- 11 Octadecane (C18)
- 12 Eicosane (C20)
- 13 Tetracosane (C24)
- 14 Octacosane (C28)
- 15 Dotriacontane (C32)
- 16 Hexatriacontane (C36)
- 17 Tetracontane (C40)
- 18 Tetratetracontane (C44)
- 19 Pentacontane (C50)
- 20 Dopentacontane (C52)
- 21 Hexacontane (C60)
- 22 Heptacontane (C70)
- 23 Octacontane (C80)
- 24 Nonacontane (C90)
- 25 Hectane (C100)

