

ASTM D7169 Crude Oil on ZB-1XT SimDist 5x0.53x0.15

Column: Zebron[™] ZB-1XT SimDist Metal, GC Cap. Column 5 m x 0.53 mm x 0.15 μ m, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 5 meters x 0.53 mm x 0.15 μ m
Order No: 7AK-G026-05
Oven Profile: 30 °C to 450 °C @ 10 °C/min for 10 min

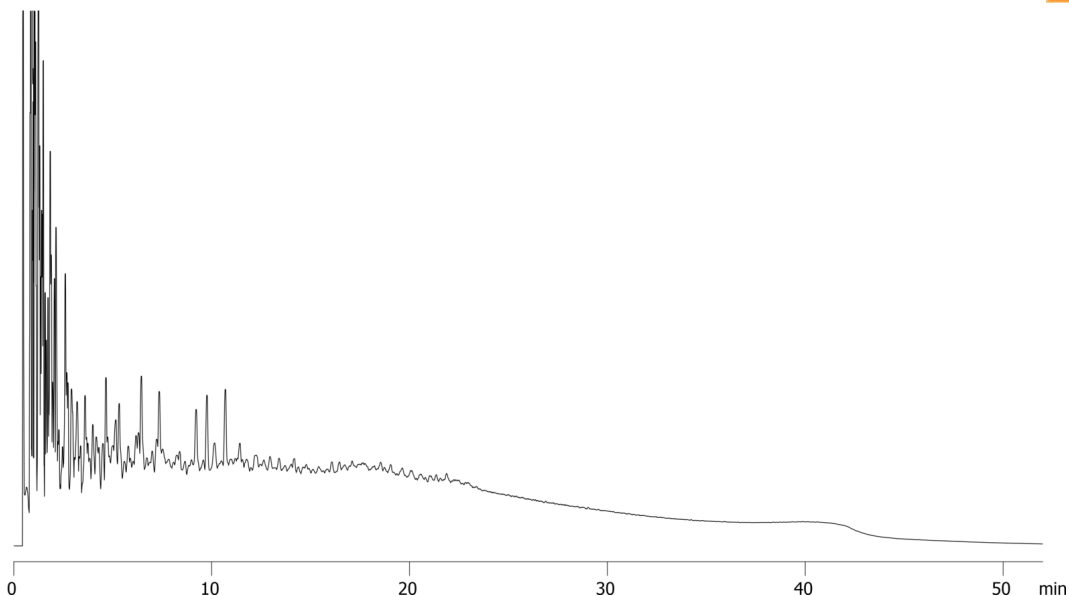
Carrier Gas: Constant Flow Helium, 7 mL/min
Injection: On-Column : 1 μ L @ 38°C
Detection: Flame Ionization (FID) (450°C)
Analyst Note: Sample was 1% in CS₂



Products used in this application:



19275

**ANALYTES:**

1 Crude oil

