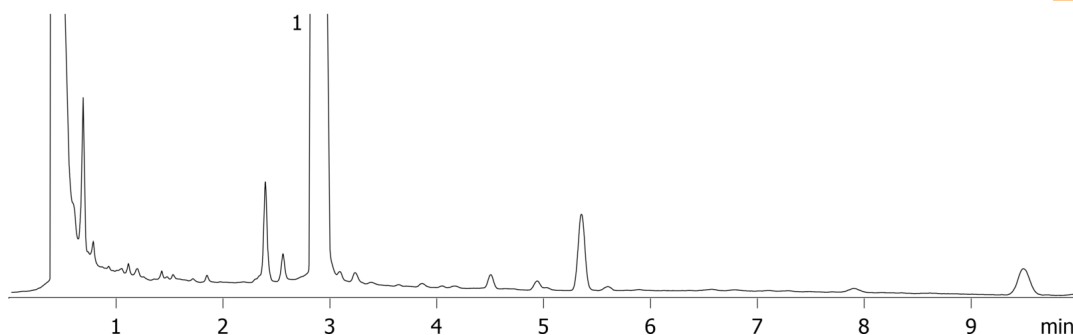


Doxylamine on a ZB-Wax plus 15x0.25x0.25 by FID**Column:** ZebronTM ZB-WAXplusTM, GC Cap. Column 15 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 15 meters x 0.25 mm x 0.25 μ m**Order No:** 7EG-G013-11**Oven Profile:** 212°C isothermal**Carrier Gas:** Constant Flow Helium, 60 cm/sec**Injection:** Split 10:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Doxylamine sample prepared as per USP 621 monograph.**Products used in this application:**

19294

**ANALYTES:****1** Doxylamine