


Local Anesthetics (Underivatized) on ZB-Drug-1, 10 m x 0.18 mm x 0.18 µm**Column:** Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:** 150 °C for 1 min to 260 °C at 25 °C/min. to 320 °C at 40 °C/min for 1 min**Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Split 10:1 1 µL @ 250°C**Detection:** Flame Ionization (FID) (320°C)**Analyst Note:** Each component is 50 ppm in methanol**Products used in this application:****ANALYTES:**

- 1 Benzocaine
- 2 Prilocaine
- 3 Lidocaine
- 4 Procaine
- 5 Tetracaine
- 6 Bupivacaine
- 7 Dibucaine

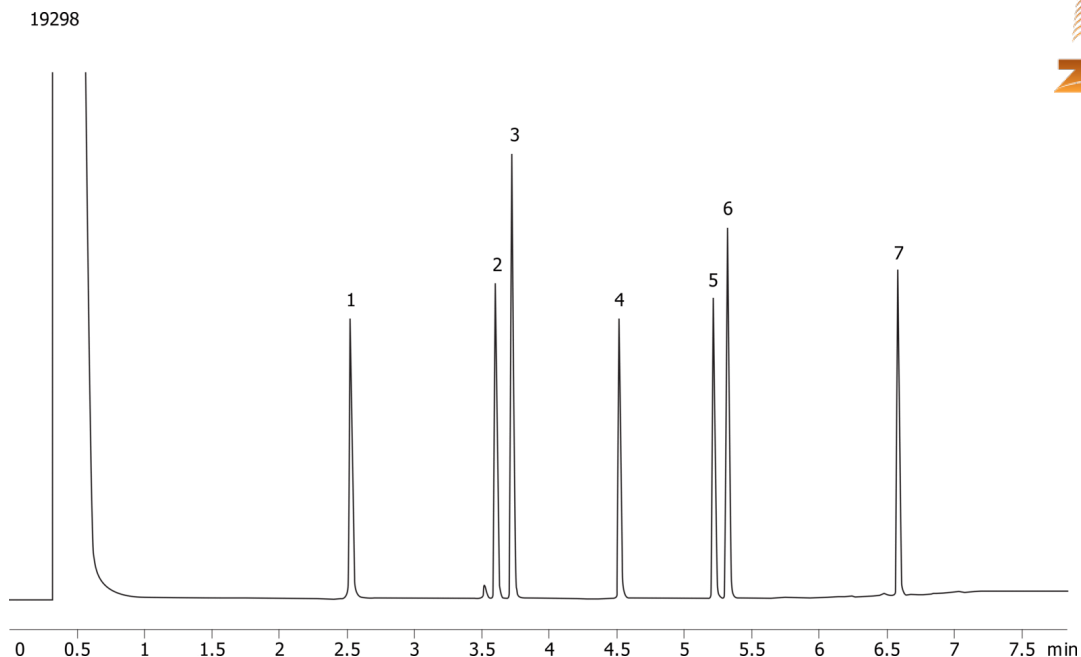