

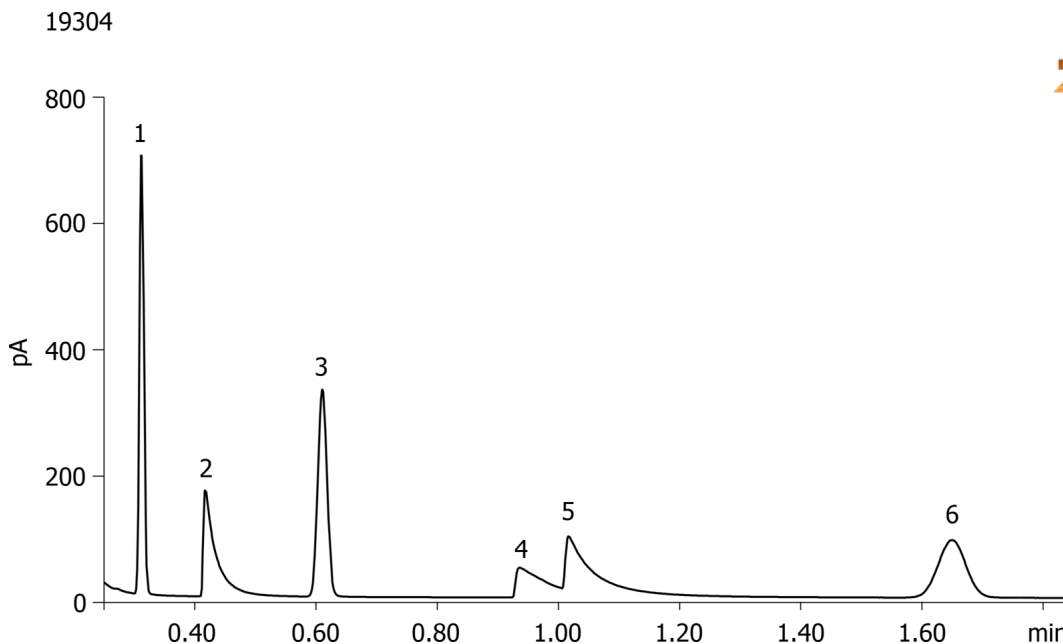
ZB-1 Test Mix on Non-Glass Infused Metal Column

Column: Zebtron™ ZB-1XT SimDist Metal, GC Cap. Column 5 m x 0.53 mm x 0.09 µm, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 5 meters x 0.53 mm x 0.09 µm
Order No: 7AK-G026-55
Oven Profile: 130 °C Isothermal

Carrier Gas: Constant Pressure Hydrogen, 1.93 mL/min
Injection: Split 100:1 1 µL @ 250°C
Detection: Flame Ionization (FID) (325°C)
Analyst Note: Sample at 250µg/mL each in Hexane



Products used in this application:

**ANALYTES:**

- 1 Undecane
- 2 4-Chlorophenol
- 3 Tridecane
- 4 1-Undecanol
- 5 Dicyclohexylamine
- 6 Pentadecane

