


ZB-1XT SimDist Test Mix**Column:** Zebtron[™] ZB-1XT SimDist Metal, GC Cap. Column 5 m x 0.53 mm x 0.15 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 5 meters x 0.53 mm x 0.15 μ m**Order No:** 7AK-G026-05**Oven Profile:** 105 °C Isothermal**Carrier Gas:** Constant Pressure Helium, 5 mL/min**Injection:** Split 20:1 0.05 μ L @ 250°C**Detection:** Flame Ionization (FID) (325°C)**Analyst Note:** Sample at 1000 μ g/mL each in Hexane**Products used in this application:**

19305

**ANALYTES:**

- 1 Undecane
- 2 Dodecane
- 3 Tridecane
- 4 1-Undecanol
- 5 Dicyclohexylamine
- 6 Pentadecane

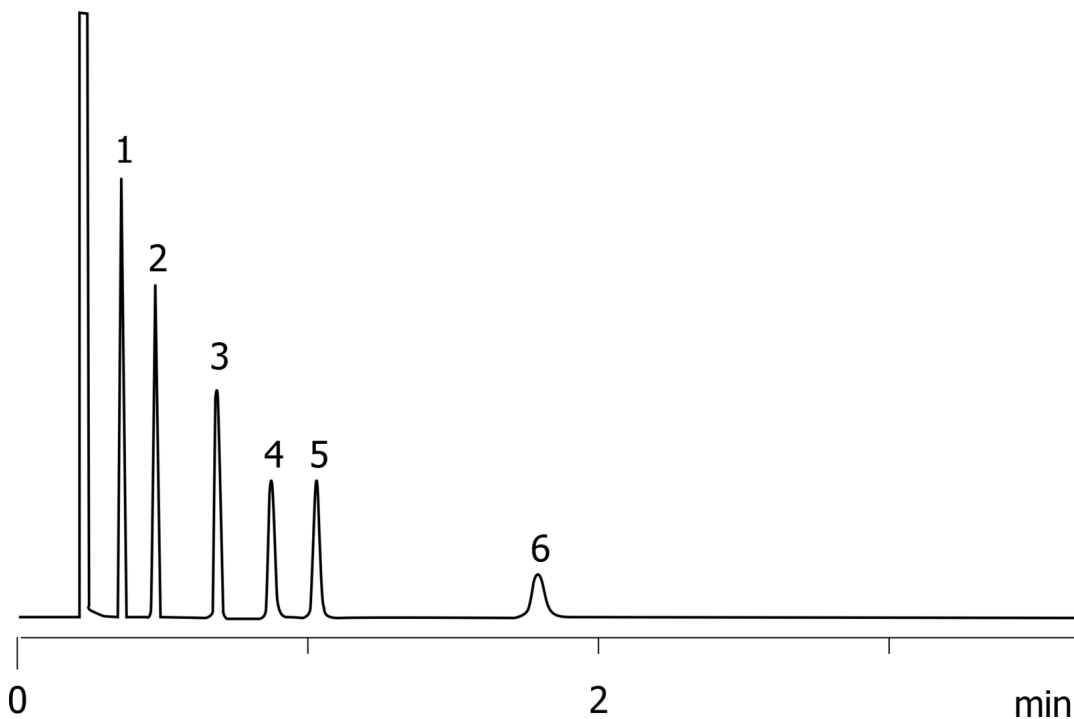